




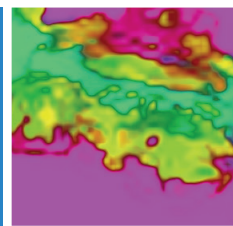

SmartSampling™: Turns Raman Imaging from Hours to Minutes

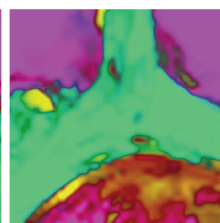

The app to speed-up
your analysis

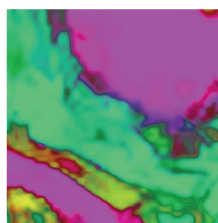

Get much more of
your sample
in less time


SmartSampling™ is HORIBA's ultimate tool for Raman Imaging. It offers incredible acquisition performance for high-resolution images, even for weak scatterers, in a fraction of the time required for standard maps.

SmartSampling™ technology uses an adaptative mapping step size approach to highlight the smallest microscopic details on the surface of your sample.

Its user friendly and intuitive approach, *SmartSampling™* supports you to make the most accurate and efficient Raman analysis. It intelligently combines precision and speed.

Use *SmartSampling™* for:

- High-resolution Raman images, with the best quality for a given time
- Large samples even with very small details

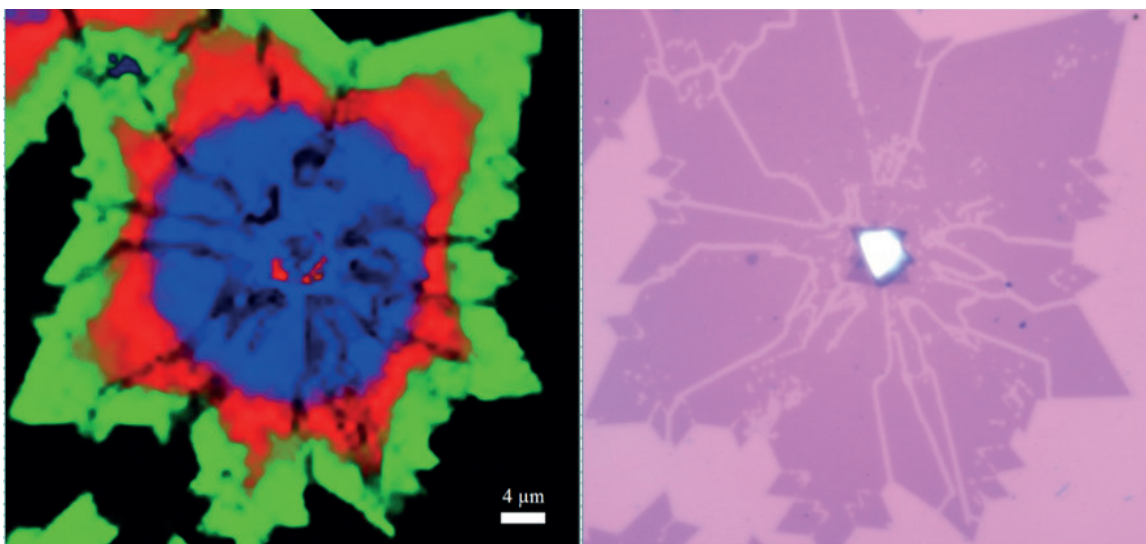


High quality Raman image of WS₂ flake: 2 layers, few layers and bulk. Total acquisition time 45 min in SmartSampling™ versus 20 h for a point by point map (sample supplier <https://2dlayer.com/>)

If an overnight point by point Raman map (100 x 100 pixels, 4 s/point) is acquired in 12 hours, an equivalent map can be obtained in a few minutes with *SmartSampling™*.


Alternatively, a 12 hour *SmartSampling™* map will allow the acquisition of a 1000 x 1000 pixels equivalent, 4 s/point map.

With a single button press, the easy to use algorithm automatically and continuously improves the quality of the complete image.

- Based on the sample image contrast *SmartSampling™* organizes the measurement points and continuously improves the image.
- The most important points are acquired first.
- The user is free to stop the scan at any time, as soon as the image quality is good enough.

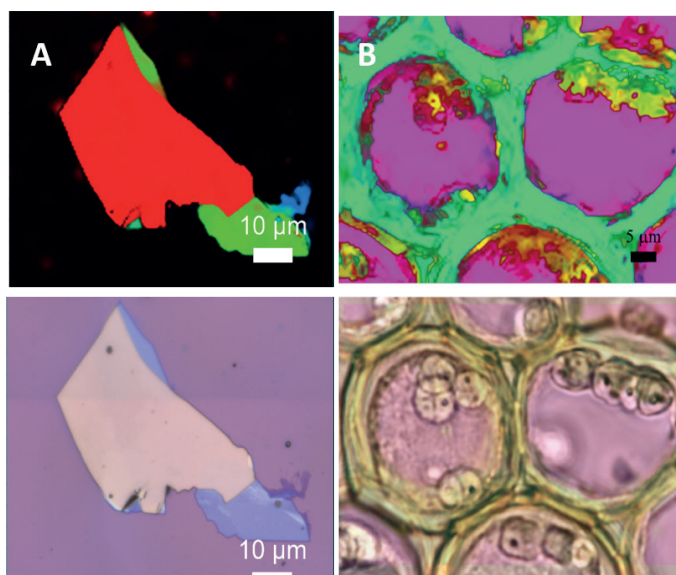

SmartSampling™ Raman Images: (A) Graphene flakes in 4 min vs 2 h 30 min in standard mode, (B) Convallaria Cells in 14 min vs 5 h in standard mode

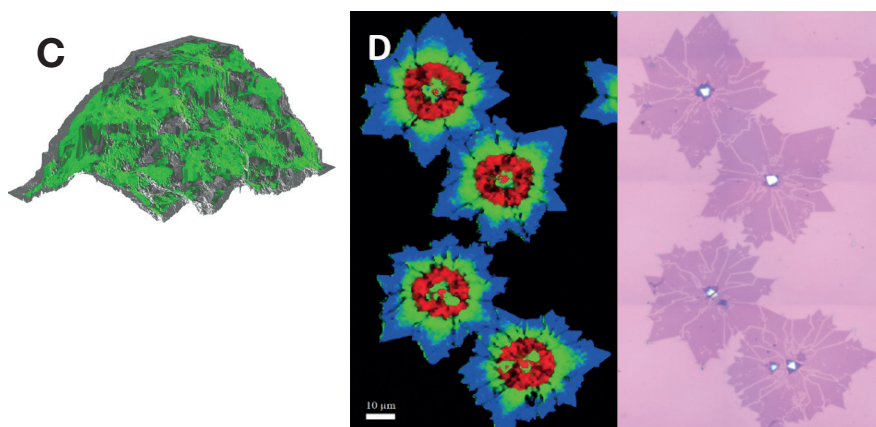

SmartSampling™ Raman Images: (C) Quartz distribution on stone with topography 3 h vs 36 h in standard mode, (D) WSe₂ flakes in 1 h 40 min vs 7 h in standard mode

HORIBA
Scientific

info.sci@horiba.com

USA: +1 732 494 8660
UK: +44 (0)1604 542 500
China: +86 (0)21 6289 6060
Taiwan: +886 3 5600606

France: +33 (0)1 69 74 72 00
Italy: +39 06 51 59 22 1
India: +91 (80) 4127 3637
Brazil: +55 (0)11 2923 5400

www.horiba.com/scientific

Germany: +49 (0) 6251 8475 0
Japan: +81(75)313-8121
Singapore: +65 (6) 745-8300
Other: +33 (0)1 69 74 72 00

20 YEARS
1919-2019
HORIBA JOBIN YVON