

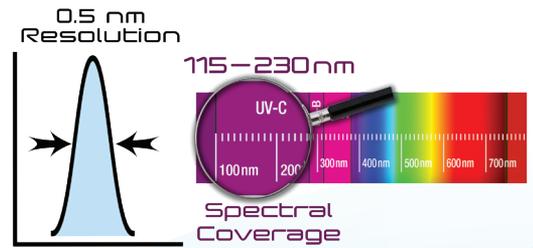
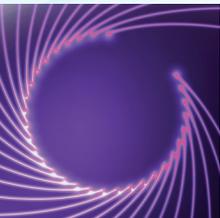


HORIBA

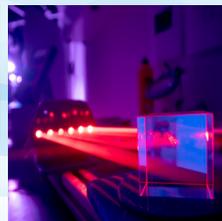

VU90 Spectrometer

Most Compact Far UV
CCD Spectrometer

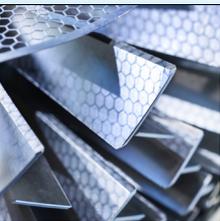

For OEM VUV
Industrial Applications

horiba.com/oem

Explore the future

HORIBA

VU90 Most Compact Vacuum UV/Far UV Back-Illuminated CCD Spectrometer

Overview

HORIBA has developed a compact size Vacuum Ultraviolet (VUV) Spectrometer with a dedicated deep-cooled CCD camera, Sincerity-VUV, all manufactured in its NJ facility. This unique system features a special CCD sensor with high sensitivity in the far UV wavelength range with capability to measure from 115 nm with a resolution better than 0.5 nm under vacuum, or high purity nitrogen. The system's internal optical cavity can be black coated for ultimate suppression of stray light.

Far UV Applications (Vacuum UV or N₂ purge)


- VUV Emission & Fluorescence
- VUV Transmission & Reflectance
- GC-VUV: Alternative to conventional gas chromatography

Examples:



- Metrology
- Environmental Monitoring
- Photolithography
- Flavors & Fragrances
- Petrology
- Food & Beverage Safety

Optical and Mechanical Layout

Based on concave holographic grating

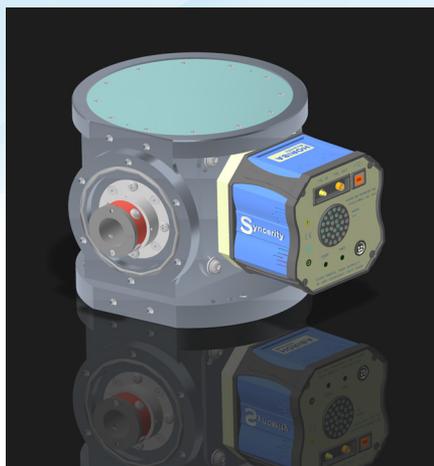

High sensitivity and very low stray light are achieved in a short focal length (90 mm) CCD spectrometer

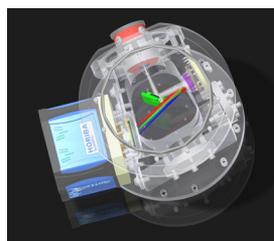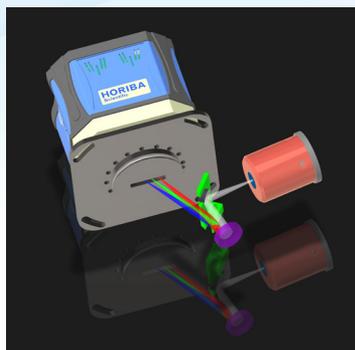

Features

First portable and most affordable VUV spectrometer

Deep-cooled CCD with high quantum efficiency

High resolution of 0.5 nm selectable resolution and throughput

Smallest VUV system footprint on the market

Ultra high throughput combined with unmatched low stray light performance

Extendable wavelength coverage from Far UV to UV-VIS-NIR


General Spectrometer Specifications

Spectral Coverage	115 - 230 nm; full coverage on CCD: 115 - 320 nm (No order sorting filter available)
Spectral Resolution	0.5 nm with 22 μm pinhole, 0.75 nm with 40 μm pinhole
Spectral Dispersion	12 nm/mm; 0.17 nm/pix
Focal Length	90 mm
VUV Optics	MgF ₂ coated grating and mirror
Vacuum or Nitrogen-purged Mode	$\sim 10^{-6}$ Torr with turbomolecular pump. Helium leak test performed at factory. N ₂ high purity nitrogen constant purge required
Options	Selection of different heights of VUV CCDs Kinematic mounts to interchange pinholes (or slits) Bare metal finish is standard; Optional super black coating reduces stray light
F/#	\sim F/10
Software	LabVIEW acquisition software for initial evaluation (DLLs provided for software integration)


Vacuum Camera Specifications (Syncerity® VUV 2048 pixel BI Sensor)*

Quantum efficiency at 20° C (See plots on next page)	Standard	Special
CCD Sensor Format	2048 x 70	Optional increased CCD sensor heights available
CCD Pixel Size	14 μm x 14 μm	
CCD Height	1 mm	
CCD QE	3 options - refer to graphs below	
Camera Window	2 options - no window for lowest VUV response or MgF ₂ window	
Deep Thermoelectric Cooling	-45° C under vacuum (10 ⁻⁶ mbar) with +25° C ambient or -25° C under N ₂ purging	
Single Pixel Well Capacity	50 000 e ⁻ /pixel (minimum); 60 000 e ⁻ /pixel (typical)	
Serial Register Full Well Capacity	250 000 e ⁻ /pixel (minimum); 500 000 e ⁻ /pixel (typical output register saturation)	
Scan Rates	45 kHz and 500 kHz	
Readout Noise (at 45 kHz and at -45° C)	10 e ⁻ (typical) to 14 e ⁻ (maximum)	
Readout Noise (at 500 kHz and at -45° C)	20 e ⁻ (typical) to 25 e ⁻ (maximum)	
Maximum Spectral Rate	20 Hz at 45 kHz scan rate; 189 Hz at 500 kHz scan rate	
Digitization	16-bit ADC	
Dynamic Range (Typical for Single Pixel)	50,000:1	
Non-linearity (Measured on Each Camera)	< 0.15% (typical) at 45 kHz (0.4% maximum)	
	< 0.20% (typical) at 500 kHz (1% maximum)	
Dark current at -45° C Note: Pixel size = 14 μm	0.1 e ⁻ /pixel/s (typical)	
Software-adjustable Gains	2, 4, and 10 e ⁻ /count at -45° C	
Communication	USB 2	
Environmental Conditions	Operating temperature 15° C to 40° C ambient	
	Relative humidity < 70% (non-condensing)	
	Storage temperature -25° C to 45° C	
Power Requirements AC/DC power supply (provided)	90-264 VAC, 47-63 Hz	

*Specifications, form factor, and spectrometer cover subject to change without notice. No LabVIEW license is needed to run our acquisition software.

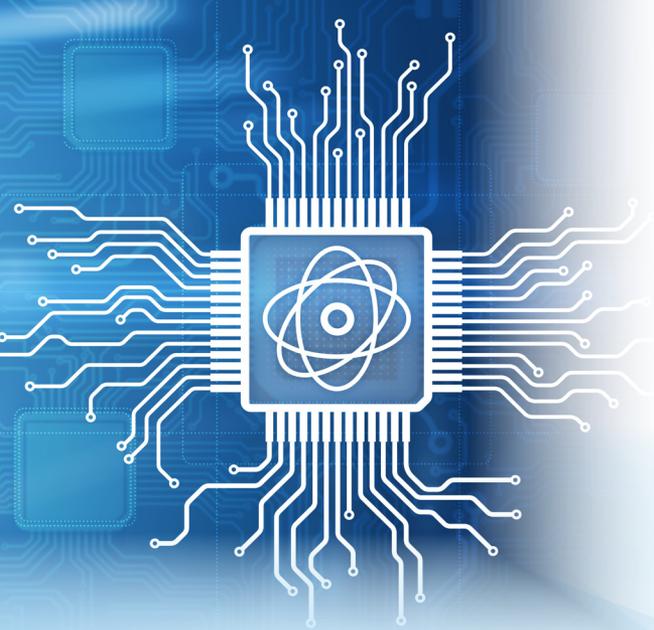


CCD Quantum Efficiency

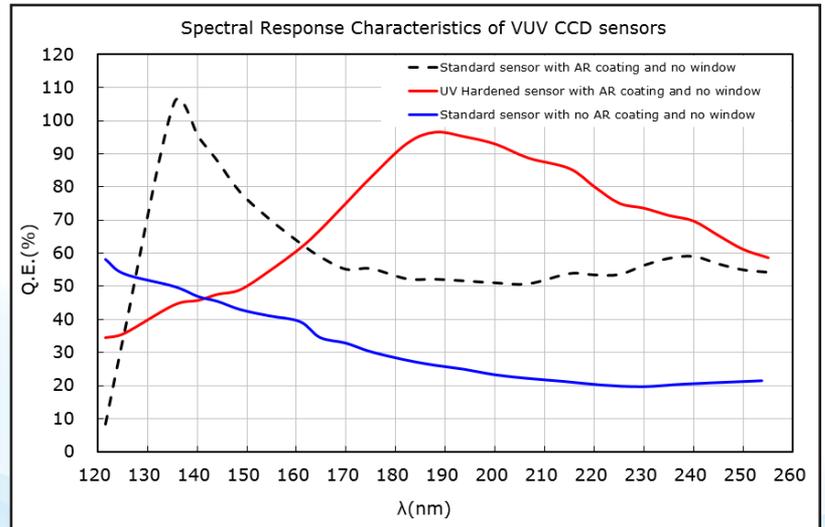


UV hardened sensor provides increased lifetime, however no warranty is provided on any sensor exposed to VUV.

The Only Commercial VUV System Using a Super Black Coating

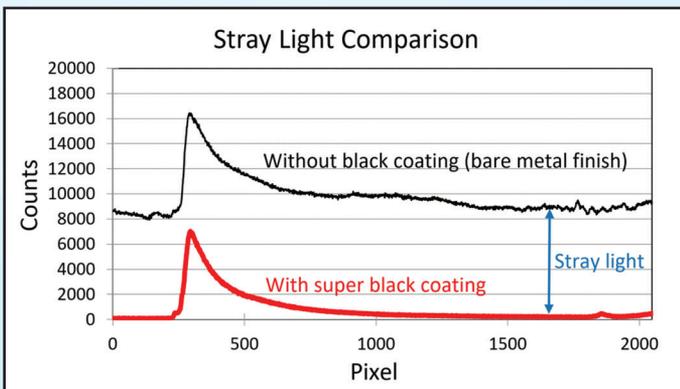


Stray light comparison between uncoated and coated cavity (Super Black) measured by far field illumination using tungsten lamp – zero order

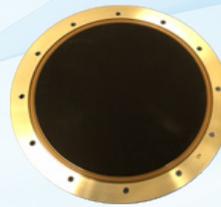


Housing Lid



Baffle treated with super black coating

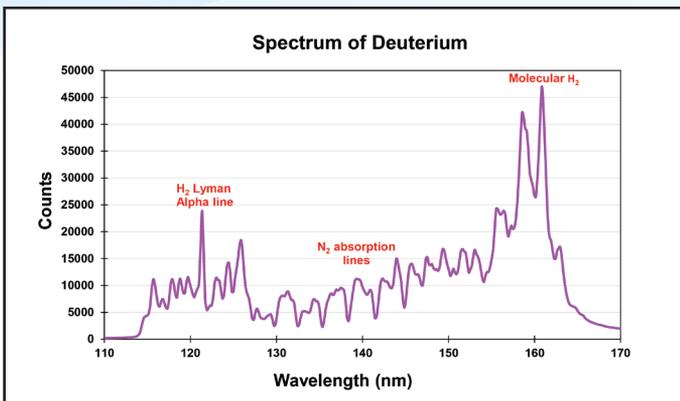



Spectrum of deuterium lamp shows enhanced sensitivity starting from 113 nm with 22 μm slit (0.45 nm resolution)

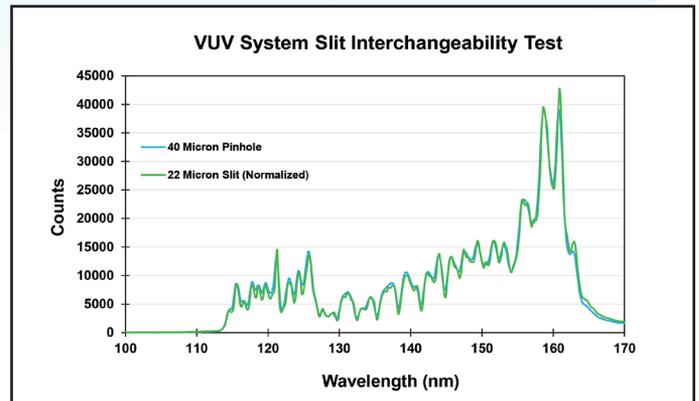



VUV system slit interchangeability test with 40 μm pinhole and 22 μm slit

Illustrations of Different Setups of FUV System

Config #1: High vacuum setup with turbo pump CCD temp: -45°C

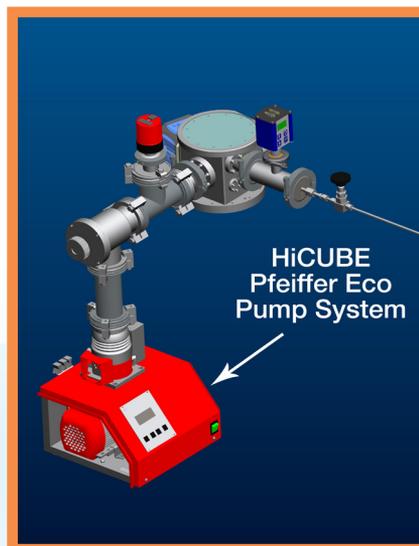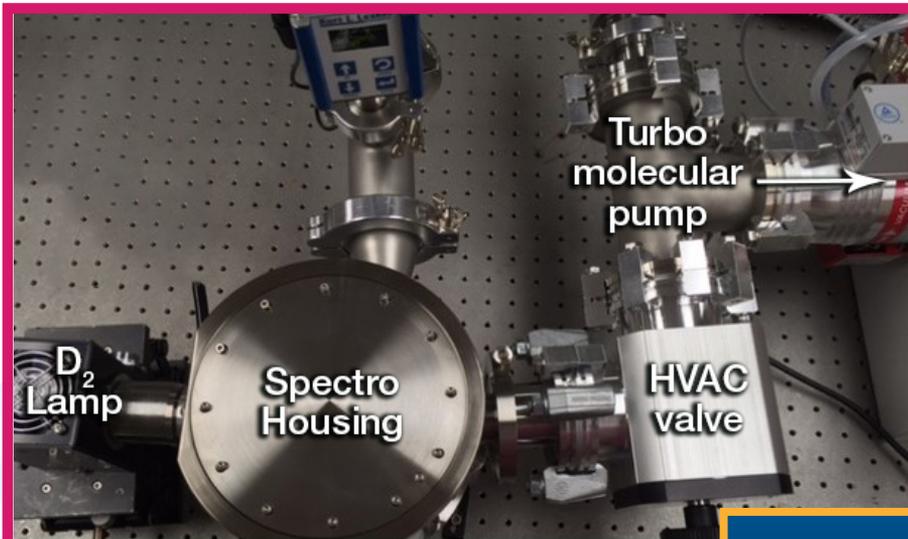

Config #2: Nitrogen-purged CCD temp: -25°C

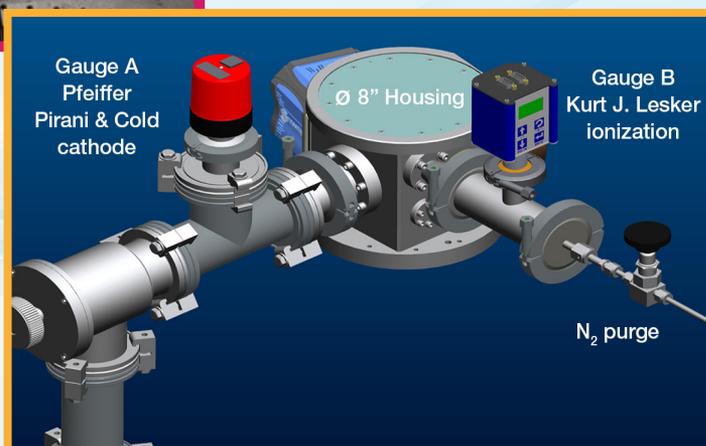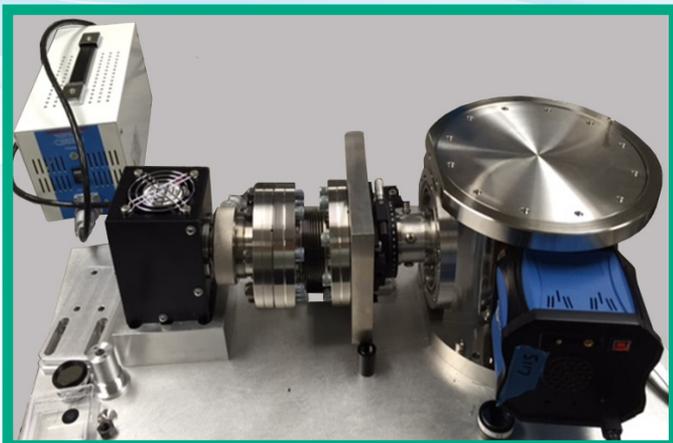


N_2 -purged set-up shown with Hamamatsu: L10366 VUV light source

Adjustable Resolution and Throughput







Resolution and throughput are adjustable by exchanging optional slit & pinhole holders.

Pinhole sizes: 22 and 40 μm diameter
 Slit sizes: 22 and 40 μm width

Kinematic interchange is achieved if all components are ordered and pre-aligned at factory.

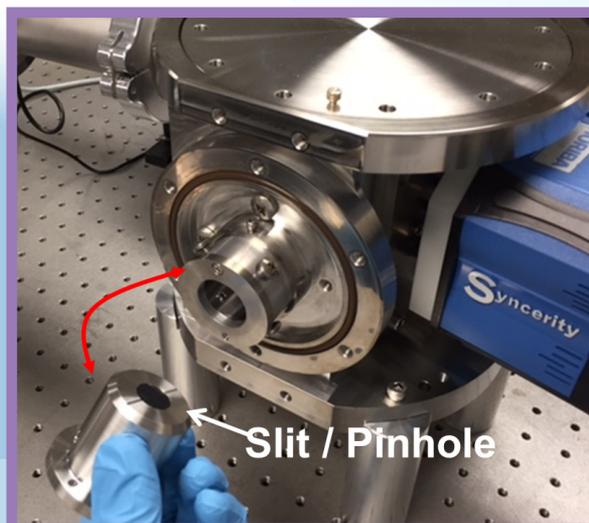

System Mechanical Drawings

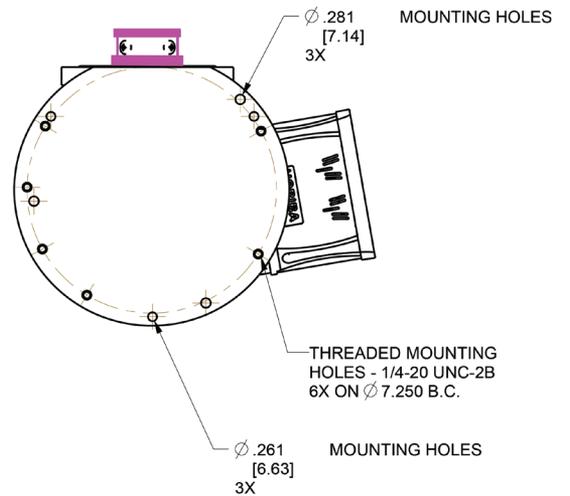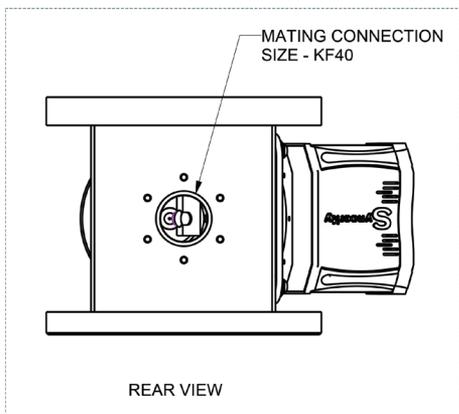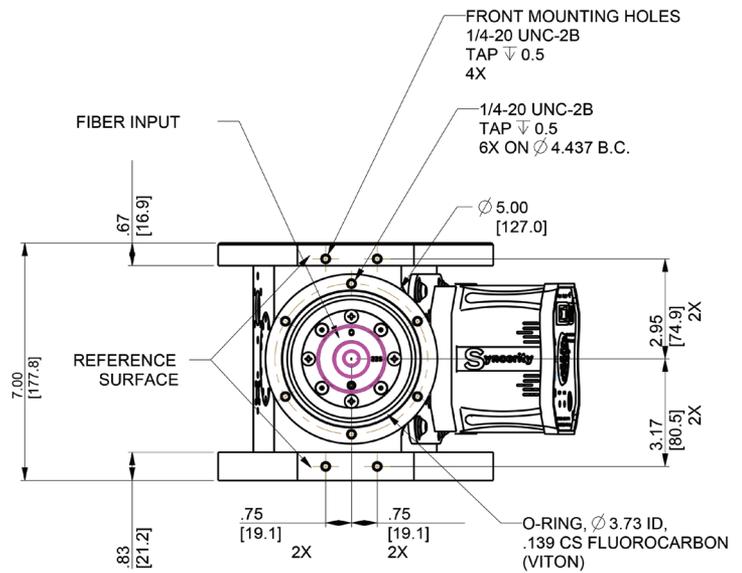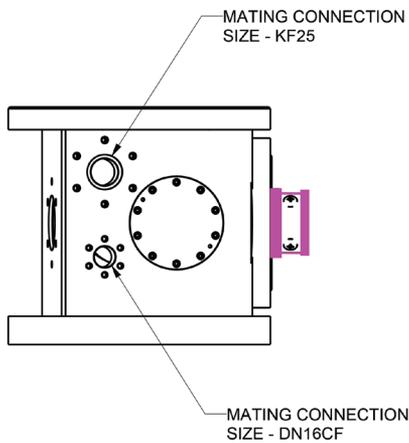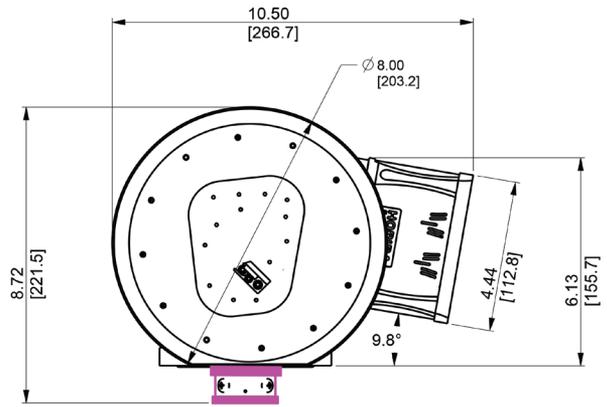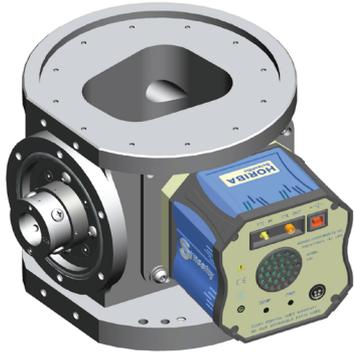


Best Selling Miniature Spectrometers for OEM Industrial Applications

Fiber-coupled USB Spectrometers:

MiniVS20 Spectrometer with Linear UV-VIS CMOS or NIR InGaAs Sensor

5 nm resolution

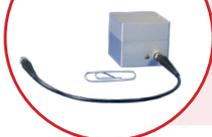

OEM hand-held spectrometer covering 190 to 1,700 nm for various low stray light applications

- Aberration-corrected concave holographic grating options
- VIS configuration featuring a 1.7" x 1.9" x 2" size combined with full F/2.3 optics for high signal-to-noise
- High throughput, compactness and long term reliability

MiniVS70 VIS Spectrometer with FI CMOS or BI CCD

1 nm resolution



NEW miniaturized VS70 configuration

- Based on high performance aberration-corrected concave gratings fitted with a custom order-sorting filter to eliminate higher orders
- Low cost combined with high performance and low stray light
- Long term opto-mechanical stability and choice of front-illuminated linear CMOS or back-illuminated CCD sensors

VS70 UV-VIS-NIR Spectrometer with Uncooled / TE-cooled CCD

1-2 nm resolution



Compact, versatile most popular VS70 OEM spectrometer and OES configurations

- Based on high performance aberration-corrected concave gratings with full F/2.3 aperture
- Affordable, high throughput, robust and stable
- Electronics drivers ranging from **USB-2 to Ethernet and EtherCAT**

CiCi-Raman-NIR with Scientific Camera Optimized for 785 nm

6 cm⁻¹ resolution



Most compact OEM Raman spectrometer with aberration-corrected holographic grating

- Covers 150-3,300 cm⁻¹
- High efficiency and low stray light
- Available in F/2.3 and in compact F/5 configurations
- -50° C deep-cooled scientific CCD camera with minimized etaloning and high NIR QE

PoliSpectra® Quad Spectrometer for Simultaneous Acquisition of 4 VIS Spectra

1 nm resolution



CCD spectrometer for simultaneous acquisition from 4 fiber inputs (470-730 nm)

- High-speed electronics (as fast as <1.5 msec readout time for 4 spectra)
- QUAD-channel high throughput system (f/2.3) and ultra-low stray light
- Industrial low-light applications from low light fluorescence to reflectance

PoliSpectra® M116 8-32 Channel MultiTrack UV-VIS-NIR CMOS Spectrometer

1 nm resolution



Fiber-coupled multi-spectra system with 8- to 32-channel simultaneous measurements

- Concentric optical design with UV extended spectral range provides minimized crosstalk
- High throughput USB-3 system featuring a fast 2D scientific BI CMOS running at 94 to 188 frames per second, acquiring 8, 16 or 32 simultaneous spectra (2048 pixels per spectrum)

PoliSpectra® 135 Multichannel or Hyperspectral Line Imager from UV to NIR


6 μm resolution





Ultra-high performance rugged spectrometer for hyperspectral imaging with a 2D sCMOS Camera

- For line-image scanning, in a push-broom hyperspectral configuration
- High throughput, USB-3 system featuring a fast 2D scientific BI CMOS with rolling shutter, running at 94 (HDR) to 188 (Standard Mode) frames per second (2048 pixels per spectrum)

OEM Spectrometer Selection Guide

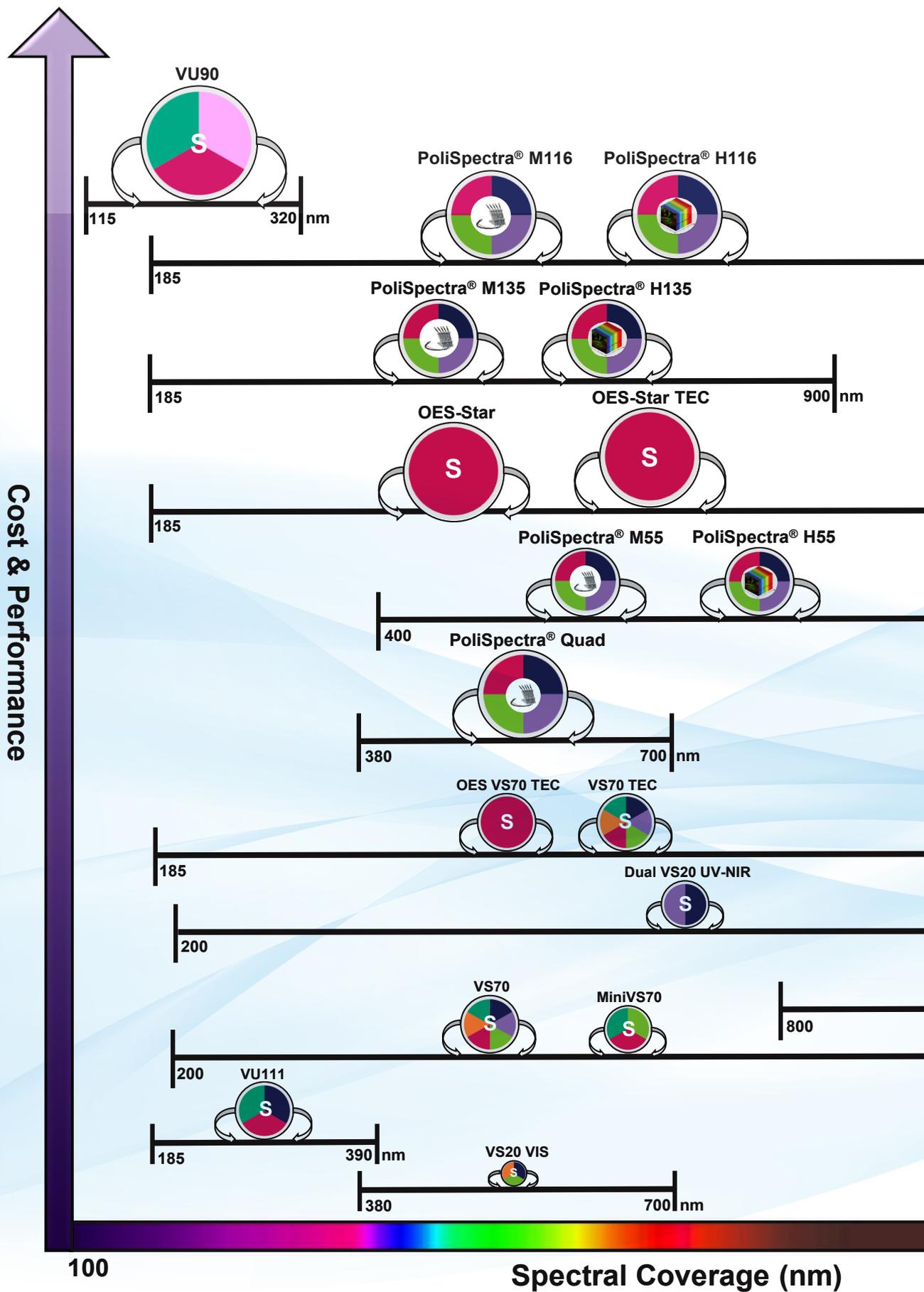

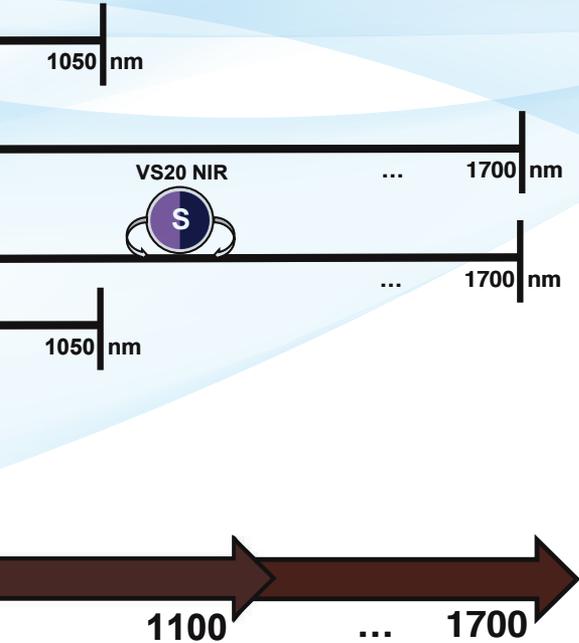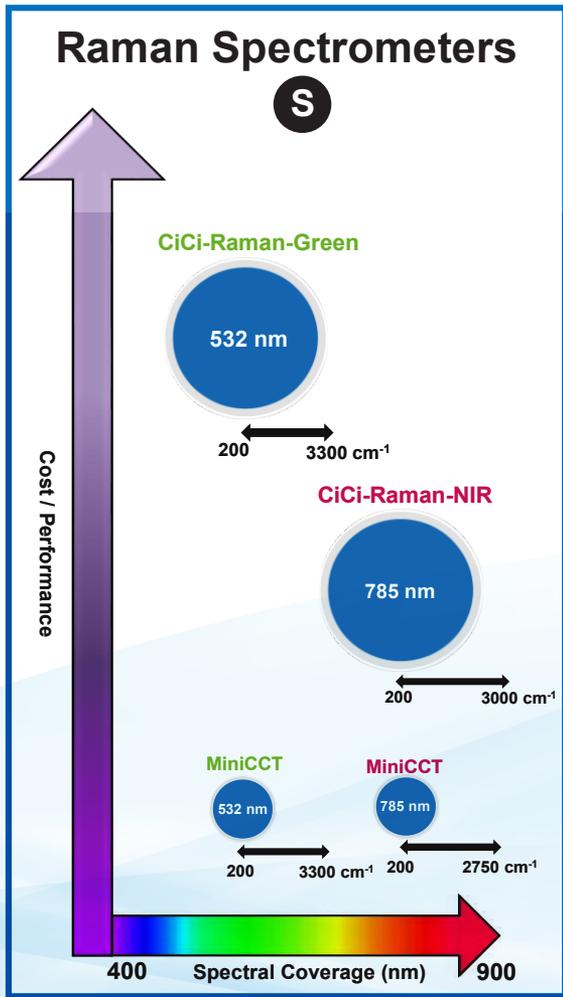

Legend

Spectrometer Input Type

S Single-fiber Input

Multi-fiber Input

Hyperspectral Imager

Technique & Application

- Emission (OES)
- Absorbance
- Fluorescence
- Reflectance
- Raman
- Metrology
- Photoluminescence
- Colorimetry

Spectrometer Size

- $\leq 10.5 \times 9 \times 7$ inches
- $\leq 8 \times 8 \times 8$ inches
- $\leq 7.5 \times 4.5 \times 5$ inches
- $\leq 7.5 \times 4.5 \times 3$ inches
- $\leq 6 \times 6 \times 7.5$ inches
- $\leq 5 \times 5 \times 5$ inches
- $\leq 5 \times 5 \times 2$ inches
- $\leq 2 \times 2 \times 2$ inches


VUV Sincerity: The latest addition to our growing family of scientific CCD Cameras with 2 available flange options

VUV Sincerity® with HORIBA VUV flange




Up to
2048 x
512
Pixels

VUV Sincerity® with Commercial flange

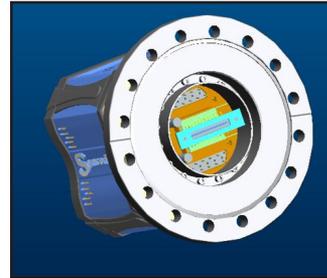

- Low Cost OEM Solution
- Available as a Stand Alone Camera
- High QE in VUV Range 58% @120 nm (no ARC)
- Attachable to High Vacuum Chamber with Viton® O-Ring, down to $<10^{-6}$ Torr
- TE-cooled to -50°C (vacuum) or -25°C with N_2 purge


Connecting to Sincerity

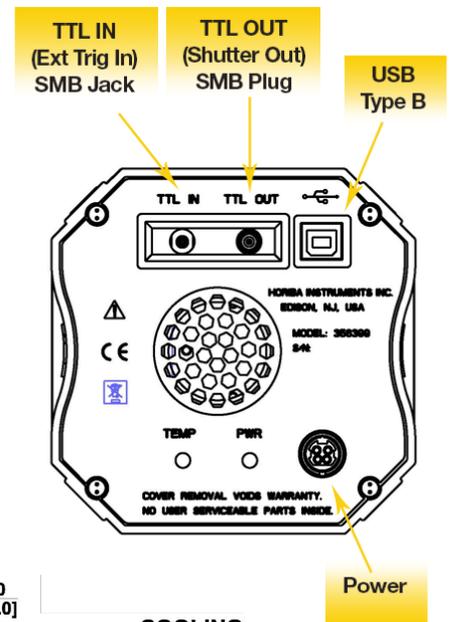

Power Interface:
Connector Type: PDP-40, Mini PWR DIN, 4-Position, STR Plug

Camera Interface:
Connector Type: USB Standard Type B

Sync I/O Connectors:
Connector Type: SMB
Input Jack: TTL IN (EXT TRIG In)
Output Plug: TTL OUT (SHUTTER Out)

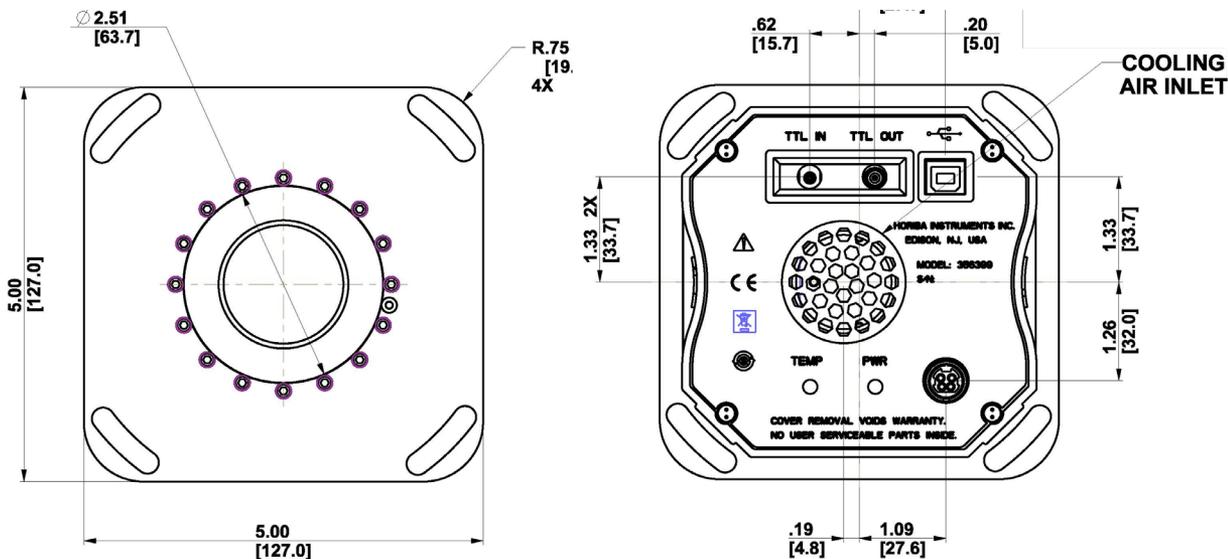

DIMENSIONS IN INCHES [MM]

OEM Philosophy and Mission

3 Centers of Excellence Dedicated to OEM Spectroscopy and Camera Solutions in US, EU, and Asia

Our mission is to provide a complete development and manufacturing experience, from optical simulations to opto-mechanical design and prototyping of spectroscopic and camera systems extending to, and including, electronics, firmware, software design and first articles.

Our products provide superior performance, reliability and stability, combined with robust cost reduction. Capable of flexible high volume production capacity in quantities of hundreds to thousands per year, we offer full confidentiality providing "Black Boxes" or private labelling, using your logo or graphics.

Unmatched customer service is provided by our exceptionally experienced workforce featuring on-time delivery and flexibility, allowing scheduling modifications.

Adhering to Copy Exactly! (CE!) processes, our fully trained staff, from engineering to manufacturing, form a dedicated OEM engineering force that supports you over the lifetime of the product.

Scientific Segment - OEM Products and Capabilities:

- Custom master optical diffraction gratings
- Diffraction grating replicas (concave, convex and flat)
- Spectrometers, optical assemblies with pre-aligned sensors (CCD, PDA, CMOS, InGaAs) using either customers' or HORIBA's OEM electronics
- OES spectrometers
- Spectroscopy systems or modular engines, such as mini fluorometers and mini Raman systems
- Single and double scanning monochromators
- Imaging spectrographs and spectrometers with CCD or CMOS cameras
- Multispectra spectrometers with multiple fiber inputs / MultiTrack spectroscopy
- Hyperspectral system with HORIBA or customer provided camera (Push-broom configurations)
- Cameras: Spectroscopic deep-cooled scientific cameras (1D and 2D CCD & InGaAs – FI and BI)
- OEM electronics for optosensors ranging from PD and PDA to CCD and CMOS sensors
- Imaging cameras: Uncooled and cooled with FI and BI high-end scientific CMOS
- VUV/FUV spectrometers and CCD vacuum and N2-purged cameras

Scientific Deep Cooled CCD, InGaAs and CMOS Cameras



Contact us in one of our centers of excellence	USA & Canada	Japan	Europe & Asia
	OEM.US@horiba.com	OEM.JAPAN@horiba.com	OEMSALES.JYFR@horiba.com
	+1 732 494 8660 Ext. 7733	+81 (75) 313 8121	+33 (0) 1 69 74 72 00