HORIBA

Instruction Manual

pH/ORP METER LAQUA-PH210 LAQUA-PH220





Preface

This manual describes the operation of the following instrument.

Brand: LAQUA

Series name: Water Quality Meter

Model: LAQUA-PH210, LAQUA-PH220

Model description: pH/ORP Meter

Be sure to read this manual before using the product to ensure proper and safe operation of the product. Also, safely store the manual so it is readily available whenever necessary. Product specifications and appearance, as well as the contents of this manual are subject to change without notice.

· Warranty and responsibility

HORIBA Advanced Techno Co., Ltd. warrants that the product shall be free from defects in material and workmanship and agrees to repair or replace free of charge, at option of HORIBA Advanced Techno Co., Ltd., any malfunctioned or damaged product attributable to responsibility of HORIBA Advanced Techno Co., Ltd. for a period of Three (3) years from the delivery unless otherwise agreed in a written statement. In any one of the following cases, none of the warranties set forth herein shall be extended:

- Any malfunction or damage attributable to improper operation
- Any malfunction attributable to repair or modification by any person not authorized by HORIBA Advanced Techno Co., Ltd.
- Any malfunction or damage attributable to the use in an environment not specified in this manual
- Any malfunction or damage attributable to violation of the instructions in this manual or operations in the manner not specified in this manual
- Any malfunction or damage attributable to any cause or causes beyond the reasonable control of HORIBA Advanced Techno Co., Ltd. such as natural disasters
- Any deterioration in appearance attributable to corrosion, rust, and so on
- · Replacement of consumables

HORIBA Advanced Techno Co., Ltd. SHALL NOT BE LIABLE FOR ANY DAMAGES RESULTING FROM ANY MALFUNCTIONS OF THE PRODUCT, ANY ERASURE OF DATA, OR ANY OTHER USES OF THE PRODUCT.

Trademarks

 Microsoft, Windows, Windows Vista are registered trademarks or trademarks of Microsoft Corporation in the United States and other countries.

Other company names and brand names are either registered trademarks or trademarks of the respective companies. (R), (TM) symbols may be omitted in this manual.

CODE:GZ0000553859B

May,2020 © 2019-2020 HORIBA Advanced Techno, Co., Ltd.

- Regulations
- EU regulations
- Conformable standards

This equipment conforms to the following standards:

 ϵ

EMC: EN61326-1

Class B, Basic electromagnetic environment

RoHS: EN50581

9. Monitoring and control instruments

Warning:

This product is not intended for use in industrial environments. In an industrial environment, electromagnetic environmental effects may cause the incorrect performance of the product in which case the user may be required to take adequate measures.

 Information on disposal of electrical and electronic equipment and disposal of batteries and accumulators

The crossed out wheeled bin symbol with underbar shown on the product or accompanying documents indicates the product requires appropriate treatment, collection and recycle for waste electrical and electronic equipment (WEEE) under the Directive 2012/19/EU, and/or waste batteries and accumulators under the Directive 2006/66/EC in the European Union. The symbol might be put with one of the chemical symbols below. In this case, it satisfies the requirements of the Directive 2006/66/EC for the object chemical. This product should not be disposed of unsorted household waste. Your correct disposal of WEEE, waste batteries and accumulators will contribute to reducing wasteful consumption of natural resources, and protecting human health and the environment from potential negative effects caused by hazardous substance in products.

Contact your supplier for information on applicable disposal methods.









Authorised representative in EU

HORIBA Europe GmbH Hans-Mess-Str.6, D-61440 Oberursel, Germany

FCC rules

FCC Compliance Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Responsible Party for FCC matter

HORIBA Instruments Incorporated Head Office 9755 Research Drive Irvine, California 92618 USA +1 949 250 4811

Note

This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at his own expense.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Korea certification

B급 기기 (가정용 방송통신기자재)

이 기기는 가정용(B 급) 전자파적합기기로서 주로 가정에서 사용하는 것을 목적으로 하며, 모든 지역에서 사용할 수 있습니다.

· Taiwan battery recycling marks



廢電池請回收

China regulation

标记的意义 Meaning of Marking

本标记适用在中华人民共和国销售电器电子产品,标记中央的数字表示环境保护使用期限的年数。(不是表示产品质量保证期间。)只要遵守这个产品有关的安全和使用注意事项,从制造日开始算起在这个年限内,不会给环境污染、人体和财产带来严重的影响。请不要随意废弃本电器电子产品。



This marking is applied to electric and electronic products sold in the People's Republic of China. The figure at the center of the marking indicates the environmental protection use period in years. (It does not indicate a product guarantee period.) It guarantees that the product will not cause environment pollution nor serious influence on human body and property within the period of the indicated years which is counted from the date of manufacture as far as the safety and usage precautions for the product are observed. Do not throw away this product without any good reason.

产品中有害物质的名称及含量

Name and amount of hazardous substance used in a product

	有害物质 Hazardous substances					
部件名称 Unit name	铅 Lead (Pb)	汞 Mer- cury (Hg)	镉 Cad- mium (Cd)	六价铬 Hexa- valent chromium (Cr (VI))	多溴联苯 Poly bromobi- phenyl (PBB)	多溴二苯醚 Poly bromo- diphenyl ether (PBDE)
本体 Main unit	×	0	0	0	0	0
电池 Battery	×	0	0	0	0	0
AC 适配器 AC adopter*1,*2	×	0	0	0	0	0
电缆 Cable ^{*2}	×	0	0	0	0	0
支架 Stand ^{*2}	0	0	0	0	0	0
打印机 Printer ^{*2}	×	0	0	0	0	0
电极 Electrode ^{*2}	×	0	×	0	0	0

本表格依据 SJ/T 11364 的规定编制。

This form is prepared in accordance with SJ/T 11364.

〇:表示该有害物质在该部件所有均质材料中的含量均在 GB/T 26572 规定的限量要求以下。

Denotes that the amount of the hazardous substance contained in all of the homogeneous materials used in the component is below the limit on the acceptable amount stipulated in the GB/T 26572.

×: 表示该有害物质至少在该部件的某一均质材料中的含量超出 GB/T 26572 规定的 限量要求。

Denotes that the amount of the hazardous substance contained in any of the homogeneous materials used in the component is above the limit on the acceptable amount stipulated in the GB/T 26572.

- *1: 本部件的环保使用期限为 10 年。 The environmental protection use period of this product is 10 years.
- *2: 选配件 Optional products

For Your Safety

For Your Safety

· Hazard classification and warning symbols

Warning messages are described in the following manner. Read the messages and follow the instructions carefully.

Hazard classification

⚠ DANGER

This indicates an imminently hazardous situation which, if not avoided, will result in death or serious injury. This is to be limited to the most extreme situations.

⚠ WARNING

This indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury.

↑ CAUTION

This indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury. It may also be used to alert against unsafe practices.

Warning symbols



Description of what should be done, or what should be followed



Description of what should never be done, or what is prohibited.

For Your Safety

Safety precautions

This section provides precautions for using the product safely and correctly and to prevent injury and damage. The terms of DANGER, WARNING, and CAUTION indicate the degree of immanency and hazardous situation. Read the precautions carefully as it contains important safety messages.

· Instrument and electrode





Do not disassemble or modify the instrument. Otherwise, it may heat up or be ignited resulting in a fire or an accident.

⚠ CAUTION

Harmful chemicals



Some electrodes are used with hazardous standard solutions. Handle them with care. The internal solution of pH electrode is highly concentrated potassium chloride (3.33 mol/L KCl). If the internal solution comes in contact with the skin, wash it off immediately. If it gets into the eyes, flush with plenty of water and then consult a doctor.



Broken glass

Broken glass may cause injury. The outer tube and tip of an electrode are made of glass. Handle them with care.



Do not use the phono jack under wet or humid conditions. Otherwise, it may cause a fire, electric shock, or breakage.

Battery

MARNING



Keep batteries out of reach of children. If someone accidentally swallows a battery, consult a doctor immediately.



If alkaline fluid from a battery gets into the eyes, do not rub the eyes, rinse with clean water immediately and then consult a doctor. Contact with alkaline fluid could cause blindness.



Do not put batteries in a fire, expose to heat, disassemble or remodel. Doing so can case fluid leakage, overheating or explosion.

Product Handling Information

Product Handling Information

Operational precautions (instrument)

- Only use the product including accessories for their intended purpose.
- Do not drop or physically impact the instrument.
- The instrument is made of solvent-resistant materials but that does not mean it is resistant to all chemicals. Do not expose the instrument in strong acid or alkali solution, or wipe with such solution.
- If the instrument is dropped into water or gets wet, wipe it using soft cloth. Do not heat to dry it.
- The instrument has a dust-proof and waterproof structure i.e., the instrument does not
 malfunction even when immersed in water of 1 m depth for 30 minutes provided that the
 electrodes are properly connected and the communication port, as well as the battery compartment, is tightly sealed. This does guarantee non-destructive, trouble-free, dust-proof,
 and waterproof performance in all situations.
- When replacing the batteries or when a serial cable connected, the instrument does not have the dust-proof and waterproof performance. The dust-proof and waterproof performance is maintained only when the covers are attached correctly.
- After replacing the batteries or removing the serial cable connected, make sure that the
 waterproof gasket attached to the cover is not deformed or discolored, or has foreign
 matter adhering to it. If the waterproof gasket is deformed, discolored or has foreign matter
 adhering to it, dust could get inside, water leaks could occur that could lead to instrument
 malfunction.
- To disconnect an electrode or serial cable, hold the connector and pull it off. If you pull at the cable, it may cause breakage.
- The phono jack communication between the instrument and a personal computer (referred to as PC in the rest of this document) may fail because of environmental conditions, such as electromagnetic noise.
- Do not replace the batteries in a dusty place or with wet hands. Dust or moisture could get inside the instrument, possibly causing instrument malfunction.
- Do not use an object with a sharp end to press the keys.
- If the power supply is interrupted while measurement data is being saved in the instrument, the data could be corrupted.
- A Ni-MH rechargeable battery can be used in this instrument.

Operational precautions (battery)

- · Do not short circuit a battery.
- Position the + and side of the battery correctly.
- When the battery has depleted or the instrument will not be used for a long time, remove the batteries.
- Of the specified battery types, make sure to use two batteries of the same type.
- Do not use a new battery together with a used battery.
- Do not use a fully charged nickel-metal hydride battery together with a partially charged battery.
- Do not attempt to charge a non-rechargeable battery.

Product Handling Information

• Environmental conditions for use and storage

- Temperature: 0 °C to 45 °C
- Humidity: under 80% relative humidity and free from condensation

• Avoid the following conditions:

- Strong vibration
- Direct sunlight
- · Corrosive gas environment
- · Locations close to an air-conditioner
- Direct wind

Transportation

When transporting the instrument, repackage it in the original package box. Otherwise, it may cause instrument damage.

Disposal

- Standard solution used for the calibration must be under neutralized before the disposal.
- When disposing of the product, follow the related laws and regulations of your country for disposal of the product.

Contents


■ Product Overview	1
Package Content	1
Key Features	2
Product components	3
■ Basic operations	7
Mode and measurement	9
■ Calibration	11
pH Calibration	11
ORP/mV Calibration	13
Temperature Calibration	15
■ Data	17
Data capture and storage	17
Data transfer	18
■ Setup	19
P1 pH setup	19
P2 Data setup	23
P3 General setup	27
P4 CLK setup	32
■ Maintenance and storage	35
 Maintenance and storage of the instrument 	35
 Maintenance and storage of pH and ORP elect 	trodes
36	
■ Error messages and trouble shooting	38
■ Appendix	
Appendix 1	
Appendix 2	
Appendix 3	47

MEMO

■ Product Overview

This section describes the package content, key features and product components of LAQUA PH200 handheld meters.

• Package Content

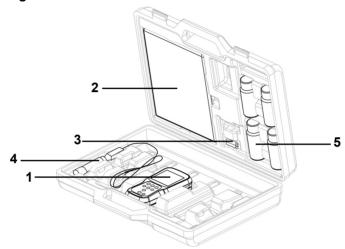

After opening the carry case, remove the meter and check for damage on the instrument and confirm that the standard accessories all exist. If damage or defects are found on the product, contact your dealer.


LAQUA PH200 Series Handheld meter and meter kit include the following items:

S.NO.	Name
1	Instrument
2	Instruction manual
3	2 AA batteries
4	Electrode
5	Calibrating solutions

Key Features


- IP67 water ingress, dust-proof, shock-resistant, anti-slip meter housing.
- Large monochrome LCD (50 x 50 mm) with white LED back lighting.
- Built-in electrode holder (up to 2 electrodes).
- · Foldable meter stand.
- · Simple user interface and single parameter display.
- 500 (for PH210) / 1000 (for PH220) data memory.
- Automatic Temperature Compensation (ATC) with temperature calibration.
- Adjustable auto shut-off time (1 to 30 minutes).
- Auto-hold / Auto stable / Real-time measurement modes with stability indicators.
- Powered by 2 x AA batteries.
- Real-time clock (only for PH220).
- PC (standard USB) / Printer (25 pin serial) connection via 2.5 mm diameter phono jack.

• Product components

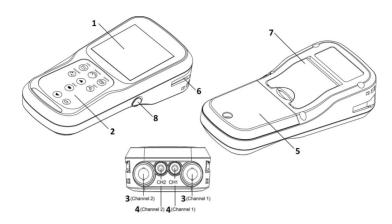

No	Name	Function
1	Monochrome LCD	Displays the measured value
2	Operation keys	Used for instrument operation
3	Electrode connector	Connect to the BNC connector of the electrode
4	Temperature connector (T)	Connect to the temperature sensor of the electrode
5	Battery cover	Open/close to insert/remove batteries
6	Electrode holder	Hold the electrode to carry with the instrument
7	Meter stand	Open stand to place the meter at an inclined position on a flat surface
8	Serial connector	Connects the serial cable and printer cable

• Display

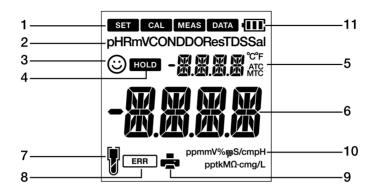

No	Name	Function
1	Status Icon	Displays the current operation mode (Setup, Calibration, Measurement and Data mode)
2	Parameters	Displays the measured parameters like pH and RmV
3	©	Stability indicator shows value is stable for documentation in Auto Stable and Auto Hold modes
4	HOLD	Appears when the measured value display is stable and fixed in auto-hold mode
5	Temperature display area	Displays the measured temperature
6	Measured value, set item display area	Displays the measured value and the set value
7	7	Indicates electrode sensitivity level
8	ERR	Indicates error situation
9	-	Indicates data being transfered to the printer or computer
10	ppmmV%mgS/cmpH pptkMΩ-cmg/L	Displays the unit for the measurement parameter
11	(III)	Displays the battery level


• Battery level display

	_
(III)	100% battery life
	50% battery life
	20% battery life
	Batteries are weak and need replacement. Refer "BATT LOW" (page 38) to solve this

• Electrode sensitivity level

¥	Electrode sensitivity above 95%(excellent)
	Electrode sensitivity between 85% to 95% (very good)
- U	Electrode sensitivity between 80% to 85%(good). Refer "SLPE ERR" (page 38) to solve this

• Keypad operation

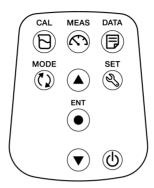

Keypad	Name	Function
9	CAL key	Switches from the measurement mode to the calibration mode. Starts calibration in the calibration mode.
5 3	MEAS key	Switches the operation mode to the measurement mode. Releases the fixed measurement value mode in the auto hold mode and begins a new measurement.
	DATA key	Switches from the measurement mode to the data mode.
(2)	MODE key	In the measurement mode, changes measurement parameters.
É	SET key	Switches from the measurement mode to the setup mode.
•	ENTER key	Determines the selection or setup. Saves data in the measurement mode and calibration mode.
•	UP key	In the setup mode, navigates between various setups. Selects preferred option in some setup screens.
•	DOWN key	Increases or decreases selected digit when entering numbers.
ტ	POWER key	Powers ON/OFF the instrument.

6

■ Basic operations

This section describes function and basic operation method of each part of LAQUA PH200 handheld meter.


• Turning on the instrument

Inserting the batteries

This instrument is operated by batteries. You can use AA alkaline batteries or AA Ni- MH chargeable batteries. Perform the following procedure to insert batteries in the instrument.

- Unscrew the battery cover on the back of the instrument counter-clock wise to unlock the battery cover.
- 2. Remove the battery cover and set the batteries inside



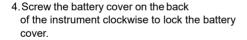




Note

- Do not replace the batteries in a dusty place or with wet hands. Dust or moisture could get inside the instrument and possibly cause an instrument malfunction.
- · Do not short-circuit a battery.
- Note polarity as shown in the battery compartment.
- When the battery has depleted or the instrument is not used for a long time, remove the batteries.
- Of the specified battery types, make sure to use two batteries of the same type.
- Do not use a new battery together with an used battery.

• Connecting an electrode

To perform calibration/ measurement, it is necessary to use the appropriate electrode for measurement parameter. Recommended electrodes for various sample are listed in our product catalog. Use the following procedure to correctly connect the electrode to the instrument:

- 1. Insert the electrode connector by fitting its groove with the connector pin of the instrument.
- 2. Turn the electrode connector clockwise by following the grooves.
- 3. Slide the connector cover on the connector.
- 4. When using a combination electrode equipped with a temperature sensor, insert the temperature jack (T) to the ATC socket on the meter.

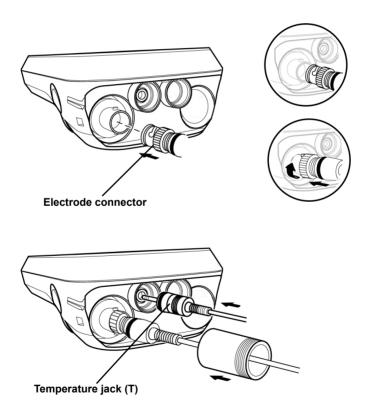

Mode and measurement


• Changing the operation mode

You can change the operation mode to four available modes depending on the purpose of use. The status icon indicates the current mode.

lcon	Name	Function
SET	Setup mode	Perform various setup functions.
CAL	Calibration mode	Performs calibration.
MEAS	Measurement mode	Performs measurement.
DATA	Data mode	Performs data setup. Displays the saved data.

You can change the operation mode using the corresponding key:

- Measurement mode: Press the 🕎 key to change to the measurement mode.
- Data mode: In the measurement mode, press the 🗒 key to change to the data mode.
- ullet Setup mode: In the measurement mode, press the $\displaystyle ullet$ key to change to the setup mode.

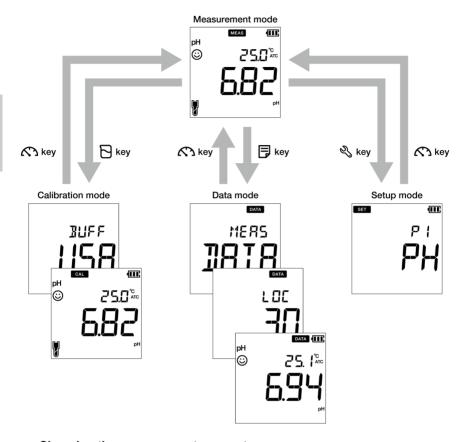

• Changing the measurement parameter

This instrument measures multiple parameters. For measurement, an electrode corresponding to the measurement parameter is required. In the measurement mode, the measurement parameter can be changed by pressing the $\ref{}$ key.



■ Calibration

This section describes the basic calibration method using LAQUA PH200 handheld meter and pH electrode.

pH Calibration

Calibration is necessary for accurate pH measurement. To perform pH calibration, follow the procedure detailed below:

Prerequisites

- Clean the pH electrode with DI (deionized) water and wipe it with tissue paper.
- Switch on the pH meter and plug in the pH electrode.
- Prepare buffer solution required for calibration.
- Keep the meter in pH measurement mode.
- Dip the pH electrode at least 3 cm in the buffer solution.

Note


- Perform two-point calibration using:
 - pH 7 and 4 for acidic sample.
 - pH 7 and 10 for alkaline sample.
- Perform three-point calibration using pH7, 4 and 10 if you are unsure of the expected sample pH value. It is recommended to calibrate with pH7 first.
- Default buffer setup is BUFF USA. If you want to change to BUFF NIST or BUFF DIN, refer to "P 1.1 Buffer selection" on page 20.



- ullet To abort an ongoing calibration process at any point of time, press the $oldsymbol{\mathcal{K}}$ key.
- It is recommended to clear the previous calibration data before performing calibration. For erasing the calibration data, refer to "P 1.3 Erase calibration data" on page 22.

Calibration

- 1. After placing the pH electrode in the buffer solution, press the \bigcap key.
- 2. The selected buffer standard appears on the meter screen and meter starts checking various calibration values with a blinking ② on screen.
- 3. Wait for the 🔾 to stabilize (stable calibration reading).
- 4. Press the **ENT** key to confirm and save calibration data.
- 5. Meter displays **DONE** indicating end of the pH calibration procedure.
- 6. Repeat for other calibration points as required.

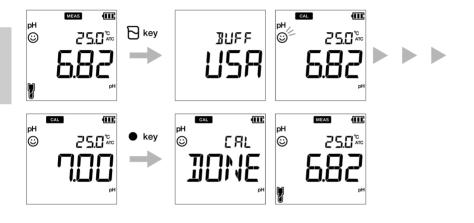

Note

If you want to know previous calibrated values, press the key when you are in the **CAL** mode. The display scrolls through the calibrated values and indicates slope and offset values.

ORP/mV Calibration

Calibration is necessary for accurate ORP measurement. To perform ORP calibration, follow the procedure detailed below:

Prerequisites

- Clean the ORP electrode with DI (deionized) water and wipe it with tissue paper.
- Switch on the pH meter and plug in the ORP electrode.
- Prepare buffer solution required for calibration.
- Ensure that the meter is in measurement mode.
- Dip the ORP electrode into the standard solution ensuring that the solution level is at least 3 cm from the electrode tip.


Note

- Absolute value measurement mode and relative value measurement mode are the two types of measurement mode available for ORP (mV) measurement.
- In absolute value measurement mode, the handheld meter displays the actual voltage value.
- In relative value measurement mode, user can adjust the absolute mV value by calibration.
 If the mV value is adjusted, the meter automatically indicates relative mV value as RmV.
 The adjustment mV is applied as an offset to the absolute mV value.
- In the relative mV mode, the absolute mV value can be adjusted by ± 200 mV.

Tip	
To abort an ongoing calibration process at any point of time, press the 🖍 key.	

Calibration

- 1. After placing the electrode in the solution, press the $\ref{2}$ key to switch to mV mode.
- 2. Press the Pkey.
- 3. Meter starts reading mV values and the 🔾 blinks until value stabilizes.
- 4. Wait for the 🛈 to stabilize (stable calibration reading).
- 5. Use the ▲ ▼ keys to adjust the mV value to your desired value.
- 6. Press the **ENT _** key to confirm and save calibration data.
- 7. Meter displays **DONE** that indicates end of the ORP/mV calibration procedure.

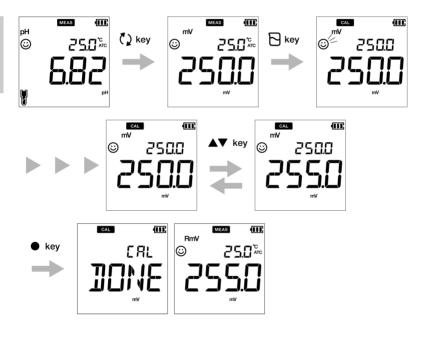

• Temperature Calibration

Temperature calibration is required to accurately match pH electrode to the meter. Check the temperature reading and if it is acceptable, no temperature calibration is required. If you need to calibrate, please follow the procedure detailed below:

Prerequisites

- Clean the pH electrode with DI (deionized) water and wipe it with tissue paper.
- Switch on the pH meter and plug in the pH electrode and temperature sensor.
- Dip the pH electrode in any buffer solution till its temperature sensor is immersed.
- Wait for 5 minutes to ensure temperature stability.


Note

- · Meter displays MTC if the temperature sensor is not plugged in and displays ATC if the temperature sensor is plugged in.
- Temperature calibration must be performed using a known temperature solution or against

a calibrated thermometer.	
Tip	
To abort an ongoing calibration process at any point of time, press the 🕰 key.	

Calibration

- 1. After placing the electrode in the solution, press the \bigcirc key.
- 2. Press the 🐧 key to switch to temperature calibration mode. Meter displays measured temperature value.
- 3. Use the **\Lambda** \textbf{V} keys to adjust the temperature to the required value.
- 4. Press the **ENT** key to save calibration data.
- 5. Meter displays **DONE** indicating end of the temperature calibration procedure.

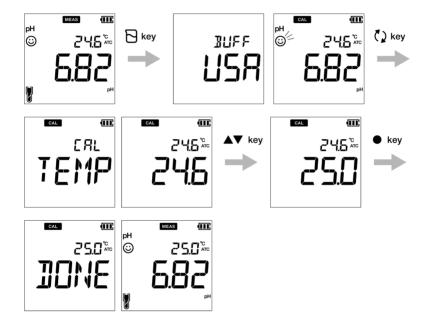

■ Data


This section describes the basic method of data storing and transferring using LAQUA PH200 handheld meters.

• Data capture and storage

In LAQUA 200 handheld meters, data measured by the instrument can be stored in the internal memory.

To save the measured data;

- Press the ENT
 key to save the displayed data.
- Meter displays the saved data for 2 seconds and then the display returns to the previous screen automatically.

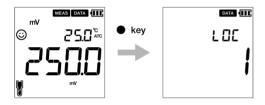

Note

- If the data storage limit reaches 500 in PH210 model or 1000 in PH220 model, memory full error occurs and MEM FULL is displayed.
- In such case, print the data or transfer necessary data to a PC (only for PH220) and delete the data from the internal memory of the instrument.

Viewing stored data

- To view stored data, press 🗒 key.
- Use ▲ ▼ keys to review different stored records.
- Press 🕎 key to return to measurement mode.



Data transfer

Transfer data to PC

Connect the instrument to a PC using the phono plug to USB cable to transfer saved data to the PC (for LAQUA PH220 only). Connect the phono jack at the instrument side to the communication port on the PC.

• Print data

To print a desired data set;

- 1. When in the measurement mode, press | key.
- 2. Use **A V** keys to view desired stored data.
- 3. Press key to print that individual data.

• Printer format- data log

Meter Model	: HORIBA PH220
Serial Number	: 123456789
SW Revision	: 1.00
User Name	:
Signature	:
Logged Data	
Location	: 2
Date	: 10 Aug 2018
Time	: 10:10:28
Mode	: pH
рН	: 7.00 pH
mV	: 0.0 mV
Temperature	: 25.0 C (MAN)
Electrode Status	: Excellent
Location	:1
Date	: 10 Aug 2018
Time	: 10:09:28
Mode	: mV
mV	: 178.0 mV
Temperature	: 25.0 C (MAN)

Tip

To print entire stored data log, refer "P 2.2 Print data setup" on page 25.

■ Setup

This section describes all the setup functions available in LAQUA PH200 handheld meters.

• P1 pH setup

Using P1 pH setup function of the meter, you can:

- · Select buffer standard
- · Set calibration alarm
- · Erase calibration data

To set the pH functions using LAQUA PH200 handheld meter, follow the procedure detailed below:

Prerequisites

Switch on the pH meter.




- Default buffer setup is BUFF USA. You can change it to BUFF NIST or BUFF DIN if required.
- Calibration alarm setup option must be used to avoid "Calibration interval alarm error" on page 38. You can set the calibration alarm for ---- day to upto 90 days, where ---indicates "no calibration alarm" has been set.
- Erasing previous calibration data is recommended for accurate calibration. Default setup is **NO** but to erase the calibration data, you need to change the setup to **YES**.

Tip	
To return to the measurement mode, press the $m{\mathcal{K}}$ key.	


P 1.1 Buffer selection

- 1. Press the 🌂 key, P1 PH screen appears.
- 2. Press the **ENT** key, **P1.1 BUFF** screen appears.
- 3. Press the **ENT** key, by default **BUFF USA** appears.
- 4. Use the ▲ ▼ keys to change the buffer standard to **BUFF NIST** or **BUFF DIN**.
- Press the ENT key, P1.1 BUFF screen appears. This indicates completion of buffer selection.

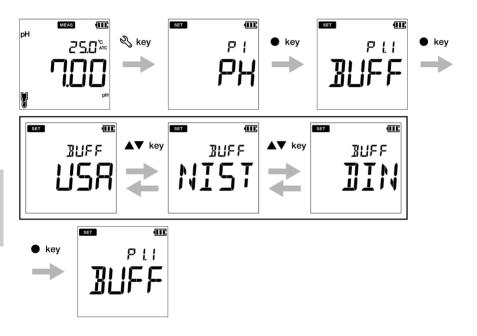


• P 1.2 Calibration alarm setup

- 1. Press the 🌂 key, P1 PH screen appears.
- 2. Press the **ENT** key, **P1.1 BUFF** screen appears.
- 3. Press the **key**, **P1.2 C.ALr** screen appears.
- 4. Press the **ENT** key, by default **DAYS** ---- appears.
- 5. Use the **A V** keys to adjust the calibration alarm interval for next calibration.
- 6. Press the **ENT •** key, **P1.2 C.ALr** screen appears. This indicates completion of calibration alarm setup.

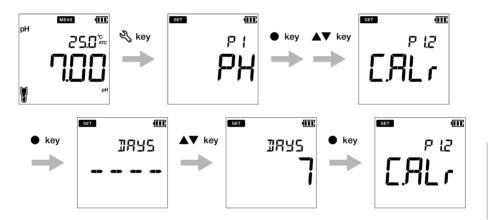

• P 1.3 Erase calibration data

- 1. Press the 🌂 key, P1 PH screen appears.
- 2. Press the **ENT** key, **P1.1 BUFF** screen appears.
- 3. Press the **k**ey, **P1.2 C.ALr** screen appears.
- 4. Press the **key**, **P1.3 C.CLr** appears.
- 5. Press the **ENT** key, **C.CLr NO** screen appears with **NO** as default setup.
- 6. Use the **A V** keys to change the setup to **YES**. This erases the calibration data.
- 7. Press the ENT key. P1.3 C.CLr screen appears. This indicates erasure of calibration data.

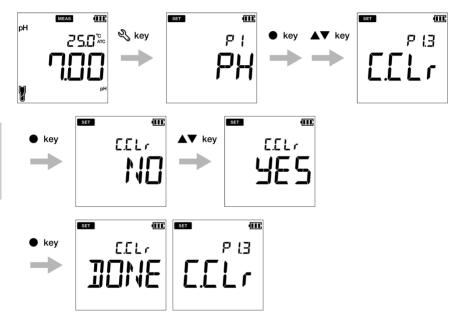

• P2 Data setup

Using P2 Data setup function of the meter, you can:

- Set data log interval
- Print data log

Prerequisites


• Erase data log

To set the data functions using LAQUA PH200 handheld meter, follow the procedure detailed below:

Switch on the pH meter.
Note
Data log interval can be set from 2 to 999 seconds.
T:-
Tip
To return to the measurement mode, press the 🔨 key.


• P 2.1 Data log interval setup

- 1. Press the \(\mathbb{N} \) key, **P1 PH** screen appears.
- 2. Press the **k**ey, **P2 DATA** screen appears.
- 3. Press the **ENT** key, **P2.1** LOG screen appears.
- 4. Press the **ENT** key, previously set log interval appears.
- 5. Use the **A V** keys to set the data log interval.
- 6. Press the **ENT** key, **P2.1 LOG** screen appears. This indicates completion of data log interval setup.

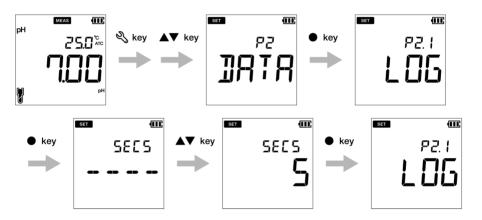


• P 2.2 Print data setup

- 1. Press the \(\mathbb{N} \) key, **P1 PH** screen appears.
- 2. Press the **k**ey, **P2 DATA** screen appears.
- 3. Press the **ENT** key, **P2.1** LOG screen appears.
- 4. Press the **key**, **P2.2 PrNT** screen appears.
- 5. Press the **ENT** key, default setup is **NO**.
- 6. Use the **A V** keys to change the setup to **YES**.
- 7. Press the **ENT** key, **P2.2 PrNT** screen appears. This indicates completion of the print data.

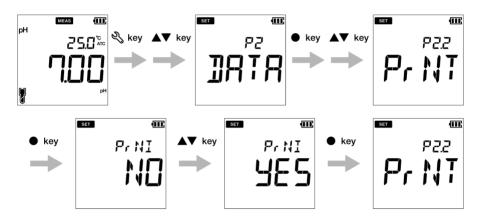

• P 2.3 Erase data

- 1. Press the 🔾 key, P1 PH screen appears.
- 2. Press **key**, **P2 DATA** screen appears.
- 3. Press the **ENT \end{align*** key, **P2.1 LOG** screen appears.
- 4. Press the **A** key, **P2.2 PRNT** screen appears.
- 5. Press the **A** key, **P2.3 D.CLR** screen appears.
- 6. Press the **ENT** key, default setup is **NO**.
- 7. Use the **\(\bigcup \)** keys set it to **YES** to erase all the data.
- 8. Press the ENT key, D.CLR DONE screen appears briefly and then P2.3 D.CLR screen appears. This indicates completion of erasure of data.

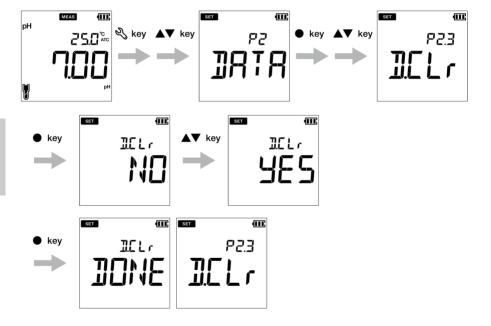

P3 General setup

Using P3 General setup function of the meter, you can:

- · Select stability mode of the meter
- · Set auto shut-off time
- Select temperature measurement
- · Reset the meter

To set the general functions using LAQUA PH200 handheld meter, follow the procedure detailed below:

Prerequisites

Switch on the pH meter.

Note

- In the calibration mode, the auto stable (AS) mode is activated. Default stability setup in measurement mode is "auto stable" (AS). If you like, you can change it to "auto hold" (AH) or "real time" (RT).
- Default auto shut-off time is 30 minutes. You can set the time from ---- to 30 minutes, where ---- indicates "no auto shut-off time" has been set and meter will be on continuously.
- Default temperature unit is °C and you can change the unit to °F.
- Default reset meter setup is NO. If you like to reset the meter, you can change it to YES.

Tip _

- Stability judgment criteria remains same for both auto stability mode and auto hold mode.
- To return to the measurement mode, press the \(\) key.


• P 3.1 Auto Stable, Auto Hold and Real Time mode setup

Auto Stable (AS) mode - the meter shows live readings ② annunciator blinks until the reading is stable.

Auto Hold (AH) mode - the meter locks the stable reading;
annunciator blinks until reading is stable and then HOLD lights up.

Real Time (RT) mode - the meter shows live readings; Both ② and HOLD annunciators are inactive.

- 1. Press the New, P1 PH screen appears.
- 2. Press **key**, **P2 DATA** screen appears.
- 3. Press **key**, **P3 GEN** screen appears.
- 4. Press the ENT key, P3.1 STBL screen appears.
- 5. Press the **ENT** key, Default the stability mode is **AS** (auto stable).
- 6. Use the \(\bigvee \) keys to change the stability mode as **AH** (auto hold) or **RT** (real time).
- Press the ENT key, P3.1 STBL screen appears. This indicates completion of the stability mode selection.

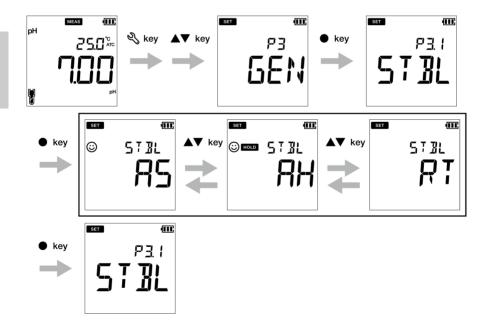


• P 3.2 Auto shut-off time setup

- 1. Press the 🌂 key, P1 PH screen appears.
- 2. Press the **A** key, **P2 DATA** screen appears.
- 3. Press the **A** key, **P3 GEN** screen appears.
- 4. Press the ENT key, P3.1 STBL screen appears.
- 5. Press the key, P3.2 A.OFF screen appears.
- 6. Press the **ENT \end{array}** key, default auto shut-off time is **30 minutes**.
- 7. Use the **A V** keys to adjust the auto off time.
- 8. Press the ENT key, P3.2 A.OFF screen appears. This indicates completion of the auto shut-off time setup.


Note

The default shut off time is 30 minutes. This can be adjusted from 1 minute to 30 minutes. If you set the display to '----' it indicates Auto Off is disabled. Meter will be on indefinitely till the user switches off the meter.

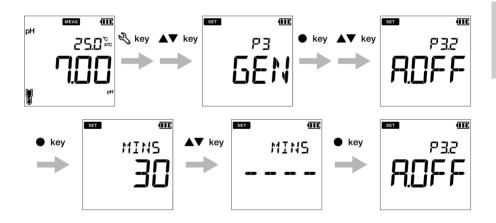


• P 3.3 Temperature unit setup

- 1. Press the key, P1 PH screen appears.
- 2. Press the **key**, **P2 DATA** screen appears.
- 3. Press the **A** key, **P3 GEN** screen appears.
- 4. Press the key, **P3.1 STBL** screen appears.
- 5. Press the key, **P3.2 A.OFF** screen appears.
- 6. Press the **key**, **P3.3°C°F** screen appears.
- 7. Press the **ENT** key, default temperature unit is °C.
- 8. Use the **A V** keys to change the unit to °F.

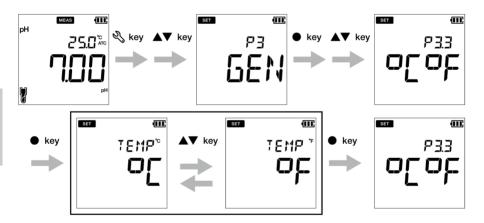

• P 3.4 Reset meter (factory default)

- 1. Press the key, P1 PH screen appears.
- 2. Press the **A** key, **P2 DATA** screen appears.
- 3. Press the **A** key, **P3 GEN** screen appears.
- 4. Press the **ENT** key, **P3.1 STBL** screen appears.
- 5. Press the **A** key, **P3.2 A.OFF** screen appears.
- 6. Press the **A** key, **P3.3** °C°F screen appears.
- 7. Press the **key**, **P3.4 r.SET** screen appears.
- 8. Press the **ENT** key, default meter re-setup is **NO**.
- 9. Use the **A V** key to set it **YES**.
- 10. Press the key. Meter displays **DONE** and automatically switches off.

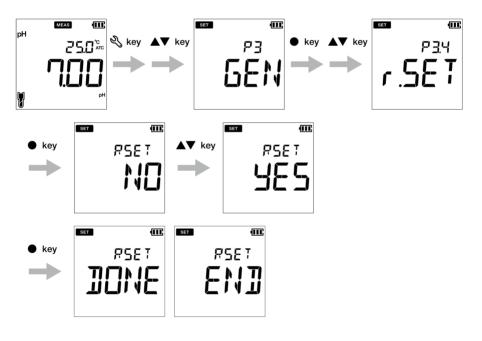

• P4 CLK setup

Real-time clock functionality is available only for LAQUA 220 meters. Using P4 Clock setup function of the meter, you can set:

- Date
- Time


To set the clock function using LAQUA PH220 handheld meter, follow the procedure detailed below:

Prerequisites

Switch on the pH meter.




- Setup date and time is necessary before using the instrument for the first time or after replacing the batteries.
- Set date and time data is captured correctly while saving data in the internal memory.

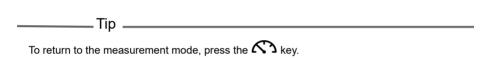


P 4.1 Date setup

- 1. Press the 🔾 key, P1 PH screen appears.
- 2. Press the **A** key, **P2 DATA** screen appears.
- 3. Press the **A** key, **P3 GEN** screen appears.
- 4. Press the **key**, **P4 CLK** screen appears.
- 5. Press the **ENT** key, **P4.1 DATE** screen appears.
- 6. Press the **ENT** key, default set year appears.
- 7. Use the \(\bigset\) keys to adjust the year.
- 8. Press the **ENT** key, default set month appears.
- 9. Use the **\(\bigvee \)** keys to adjust the month.
- 10. Press the **ENT \end{array}** key, default set day appears.
- 11. Use the A wkeys to adjust the day.
- 12. Press the **ENT** key, **P4.1 DATE** screen appears. This indicates completion of the date setup.

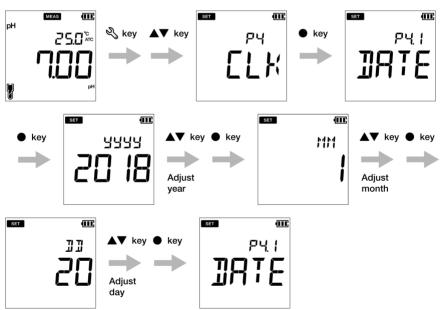

• P 4.2 Time setup

- 1. Press the 🔾 key to switch to the setup mode, P1 PH screen appears.
- 2. Press the **k**ey, **P2 DATA** screen appears.
- 3. Press the **key**, **P3 GEN** screen appears.
- 4. Press the **\(\)** key, **P4 CLK** screen appears.
- 5. Press the ENT key, P4.1 DATE screen appears.
- 6. Press the **key**, **P4.2 TIME** screen appears.
- 7. Press the **ENT** key, default set hour appears.
- 8. Use the **\(\bigvee \)** keys to adjust the hour.
- 9. Press **ENT** key, default set minute appears.
- 10. Use the **\(\bigvere** keys to adjust the minute.
- 11. Press ENT key, P4.2 TIME screen appears. This indicates completion of the time setup.

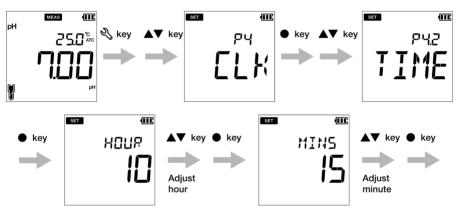



■ Maintenance and storage

This section describes maintenance of LAQUA PH200 handheld meters, pH and ORP electrodes used with the meter.

Maintenance contract

Please contact your dealer for the product maintenance contract.

Maintenance and storage of the instrument

How to clean the instrument

- If the instrument is dirty, wipe it gently with a soft dry cloth. If it is difficult to remove the dirt, wipe it gently with a cloth moistened with alcohol.
- •The instrument is made of solvent resistant materials but is not resistant to all chemicals.Do not dip the instrument in strong acid or alkali solution, or wipe it with such solutions.
- Do not wipe the instrument with polishing powder or other abrasive compound.

• Environmental conditions for storage

- Temperature: 0 °C to 45 °C
- Humidity: under 80% relative humidity and free from condensation

• Avoid the following conditions:

- · Dusty place
- Strong vibration
- · Direct sunlight
- · Corrosive gas environment
- · Close to an air-conditioner
- Direct wind

• Maintenance and storage of pH and ORP electrodes

This section describes an overview of the procedures for maintenance and storage of pH electrode.

How to clean the electrodes

When the tip of an electrode (responsive membrane and liquid junction) becomes dirty, the response time may slow or an error may occur in the measurement results. To avoid such error, clean the electrode. For dirt that cannot be washed off by pure water (or deionized water), use the cleaning solution indicated below depending on the type of dirt. After cleaning, rinse the electrode with pure water (or deionized water).

However for pH and ORP electrodes, different cleaning solutions should be used to clean different types of dirt.

For pH electrode

Type of dirt	Cleaning solution	
General	Diluted neutral cleaning solution	
Oil	Alcohol, or diluted neutral cleaning solution	
Inorganic substance	1 mol/L HCl or electrode cleaning solution	
Protein	Cleaning solution including protein-removing enzyme	
Alkali	Dip in 1 mol/L HCl or electrode cleaning solution for 1h to 2 h	

For ORP electrode

Type of dirt	Cleaning solution	
General	Diluted neutral cleaning solution (General dish washing liquid works reasonably well.)	
Oil		
Inorganic substance	Immerse dilute nitric acid (1:1 nitric acid)	

• Daily storage of the electrodes

If the electrode becomes dry, the response will be slow. Store in a moist atmosphere. Follow the steps below to properly store the electrode:

- 1. Wash the electrode well with pure water (or deionized water) to remove sample completely, and close the internal solution filler port.
- 2. Wash the inside of the protective cap with pure water (or deionized water), then add enough pure water (or deionized water) to soak the sponge.
- 3. Attach the protective cap.



When the electrode will not be used for a long period, store it by following the electrode storage procedure detailed above.

■ Error messages and trouble shooting

• Error message

This section describes the causes of typical errors and the actions to be taken to resolve respective errors.

If ERR is displayed while you are using the instrument, check the error, its cause and action to be taken in the error list below:

Meter display	ERR description	Cause of error and How to solve the problem
BATT LOW	Low battery	Battery power is low. Please replace with new batteries.
OFFS ERR	Offset voltage error	Electrode is dirty or reference junction is clogged. Clean the electrode.
SLPE ERR	Slope error	Electrode sensitivity is low. Please clean and recalibrate with fresh standard solution. If the problem persists, replace the electrode with new one.
BUFF ERR	Can not auto recognize standard solution	The instrument cannot identify the standard solution. Check the calibration solution and use fresh one if required.
	Calibration interval alarm error	Exceeds the calibration interval setup. Calibrate the meter.
MEM FULL	Memory data full	The number of the data saved has exceeded the specified number of items. Print or transfer the data. Or, clear stored data.
ERR	If user selects the enter key before stable in calibration mode	 ▶ key is pressed before the calibration value has stabilized. Wait for the value to be stable and then press the ▶ key.

• Trouble shooting

This section describes causes and actions to take for problems that customers frequently ask.

The indicated value fluctuates

< Problem with the electrode >

Cause	How to solve problem
The electrode is dirty.	Clean the electrode.
The electrode is broken.	Replace the electrode.
The level of reference electrolyte gel.	Replace the electrode.

< Problem with the instrument >

Cause	How to solve problem
There is a motor or other device causing electrical interference.	Measure at a place where no influence from induction is given. Ground all AC-powered equipment.
The electrode is not connected correctly.	Connect the electrode properly.

< Problem with the sample >

Cause	How to solve problem
Electrode is not immersed enough to cover liquid junction.	The electrode must be immersed up to the liquid junction. As a guide, immerse to at least 3 cm from the tip of the electrode.
The stability of electrode is affected by the sample solution.	It is important to select an electrode that is appropriate for the sample. Consult your dealer. To confirm an electrode that is appropriate for the sample, check the pH electrode selection guide in our catalogue, or refer to our website.

The response is slow

Cause	How to solve problem
The electrode is dirty.	Clean the electrode.
The electrode is broken.	Replace the electrode.
The response of electrode is affected by the sample solution.	It is important to select an electrode that is appropriate for the sample. Consult your dealer. To confirm an electrode that is appropriate for the sample, check the pH electrode selection guide in our catalogue, or refer to our website.

The indicated value does not change

Cause	How to solve problem
The electrode is broken.	Replace the electrode.
The electrode is not connected correctly.	Connect the electrode correctly.
Keys are locked.	Turn OFF the power, remove the batteries, and then turn ON the power again.
The instrument is in HOLD state.	Push the MEAS, and cancel the HOLD state.
Instrument defect	Consult your dealer.

The measured value is out of the measurement range

When the measured value is below the display range, "Ur" appears. When the measured value is over the display range, "Or" appears.

Cause	How to solve problem
Sample is out of the measurement range.	Use a sample within the measurement range.
Electrode is not immersed enough to cover liquid junction.	The electrode must be immersed up to the liquid junction. As a guide, immerse to at least 3 cm from the tip of the electrode.
The electrode cable is broken.	Replace the electrode.
Calibration is not performed or performed incorrectly.	Perform calibration correctly.
Instrument defect	Check as explained below.

• How to check for instrument defect

Short the metal part of the outer tube to the center pin of the electrode connector of the corresponding channel of the instrument. If "Ur" or "Or" appears in this condition, consult your dealer.



Repeatability of the measured value is poor.

Cause	How to solve problem
Effect of the sample solution.	Repeatability becomes poor when the pH of the sample changes over time.
The electrode is dirty.	Clean the electrode.
The electrode is broken.	Replace the electrode.
The internal solution of the electrode is partially depleted or contaminated.	Replace the electrode.

Nothing appears when the power is turned ON

Cause	How to solve problem
Power is not supplied.	Insert batteries.
Battery polarity (+, -) is reversed.	Insert the batteries with the polarity (+, –) correctly oriented.
Battery life is low.	Replace the batteries.
Instrument defect.	Consult your dealer.

Swelling of keypad

Cause	How to solve problem
Using the instrument at high elevation or other location where the air pressure is different from sea level.	To eliminate the pressure difference between the inside and outside of the instrument, briefly open and then close the serial connector and battery cover. After opening, correctly close the cover to maintain dust and water proofing.
Instrument defect	Consult your dealer.

Part of the display is missing

Cause	How to solve problem	
Instrument defect.	Check the display by switching ON the instrument when all the LCD segments are lit.	

■ Appendix

Appendix 1

This section describes technical information and option for LAQUA PH200 handheld meters.

pH calibration can be performed according to several buffer standards. The most common standard is the US buffer standard. The default setup for the LAQUA PH200 Series meters is US buffer standard. Alternative standards that can be chosen with the LAQUA 200 meter are NIST and DIN.

The pH buffers are temperature dependent i.e. the pH value changes with change in temperature. The LAQUA PH200 is intelligent to detect the temperature and pH value associated with buffer when calibration is performed. The meter remembers all the temperature vs. pH value for all the standards.

The pH vs. temperature values for the various standards are listed below:

<USA>

Temp. (°C)	pH 1.68	pH 4.01	pH 7.00	pH 10.01	pH 12.46
0	1.67	4.01	7.12	10.32	
5	1.67	4.01	7.09	10.25	13.25
10	1.67	4.00	7.06	10.18	13.03
15	1.67	4.00	7.04	10.12	
20	1.68	4.00	7.02	10.06	12.64
25	1.68	4.01	7.00	10.01	12.46
30	1.69	4.01	6.98	9.97	12.29
35	1.69	4.02	6.98	9.93	
40	1.70	4.03	6.97	9.89	11.99
45	1.70	4.04	6.97	9.86	
50	1.71	4.06	6.97	9.83	11.73
55	1.72	4.08	6.97	9.81	

< NIST >

Temp. (°C)	pH 1.68	pH 4.01	pH 6.86	pH 9.18	pH 12.46
0	1.67	4.00	6.98	9.46	
5	1.67	4.00	6.95	9.39	13.25
10	1.67	4.00	6.92	9.33	13.03
15	1.67	4.00	6.90	9.27	
20	1.68	4.00	6.88	9.22	12.64
25	1.68	4.01	6.86	9.18	12.46
30	1.69	4.01	6.85	9.14	12.29
35	1.69	4.02	6.84	9.10	
40	1.70	4.03	6.84	9.04	11.99
45	1.70	4.04	6.83	9.04	
50	1.71	4.06	6.83	9.01	11.73
55	1.72	4.08	6.83	8.99	

< DIN >

Temp. (°C)	pH 1.09	pH 3.06	pH 4.65	pH 6.79	pH 9.23	pH 12.75
0	1.08	3.10	4.67	6.89	9.48	13.37
5	1.09	3.10	4.66	6.87	9.43	13.37
10	1.09	3.10	4.66	6.84	9.37	13.37
15	1.09	3.08	4.65	6.82	9.32	13.17
20	1.09	3.07	4.65	6.80	9.27	12.96
25	1.09	3.06	4.65	6.79	9.23	12.75
30	1.10	3.05	4.65	6.78	9.18	12.61
35	1.10	3.04	4.65	6.77	9.13	12.45
40	1.10	3.04	4.66	6.76	9.09	12.29
45	1.11	3.04	4.67	6.76	9.04	12.14
50	1.11	3.04	4.68	6.76	9.00	11.98
55	1.11	3.04	4.69	6.76	8.96	11.84

Note

Calibration is performed using Nernst's equation with the above values.

• Appendix 2

Printer format - Measurement

рΗ

Meter Model	: HORIBA PH220
Serial Number	: 123456789
SW Revision	: 1.00
Date	: 20 Aug 2018
Time	: 10:10:28
Mode	: pH
рН	: 7.00 pH
mV	: 0.0 mV
Temperature	: 25.0 C (MAN)
Electrode Status	: Excellent
User Name	:
Signature	:

mV

Meter Model	: HORIBA PH220
Serial Number	: 123456789
SW Revision	: 1.00
Date	: 20 Aug 2018
Time	: 10:10:28
Mode	: mV
mV	: 0.0 mV
Temperature	: 25.0 C (MAN)
User Name	:
Signature	:

Relative mV

Meter Model	: HORIBA PH220
Serial Number	: 123456789
SW Revision	: 1.00
Date	: 20 Aug 2018
Time	: 10:10:28
Mode	: R. mV
R. mV	: 3.0 mV
Offset	: -3.0 mV
Temperature	: 25.0 C (MAN)
User Name	:
Signature	:

Printer format - Data log

Meter Model	: HORIBA PH220
Serial Number	: 123456789
SW Revision	: 1.00
User Name	:
Signature	:
Logged Data	
Location	: 2
Date	: 10 Aug 2018
Time	: 10:10:28
Mode	: pH
рН	: 7.00 pH
mV	: 0.0 mV
Temperature	: 25.0 C (MAN)
Electrode Status	: Excellent
Location	: 1
Date	: 10 Aug 2018
Time	: 10:09:28
Mode	: mV
mV	: 178.0 mV
Temperature	: 25.0 C (MAN)

Printer format - Calibration

PΗ

Meter Model	: HORIBA PH220
Serial Number	: 123456789
SW Revision	: 1.00
Date	: 20 Aug 2018
Time	: 10:10:28
Calibration data	
Cal Date	: 20 Jun 2018
Cal Time	: 10:10:10
Cal Points	: 4.01, 7.00, 10.01
Offset	: 0.0 mV
Avg Slope	: 98.2 %
Cal Temp.	: 25.0 C (ATC)
Electrode Status	: Excellent
User Name	:
Signature	:

• Appendix 3

Model	PH210	PH220	
Model	pH/ORP/Temp (°C/°F)		
pH Range	-2.00 to 16.00 pH		
Resolution	0.01 pH		
Accuracy	±0.0	1 pH	
Calibration Points	USA & NIST (Up t	o 5), DIN (Up to 6)	
pH Buffer Groups	USA, N	IST, DIN	
ORP Range	±200	0 mV	
Resolution	0.1 mV (< ±1000 mV), 1 mV (≥ ±1000mV)	
Accuracy	±0.3 mV (< ±1000 mV), 0.3	3% of reading (≥ ±1000mV)	
Calibration Option	Ye	es	
Temperature Range	-30.0 to 130.0 °C	/ -22.0 to 266.0 °F	
Resolution	0.1 °C	C / °F	
Accuracy	± 0.5 °C	/ ± 0.9 °F	
Calibration Option	Yes		
Memory	500	1000	
Auto Data Log	•	•	
Real-time Clock	-	•	
Date & Time Stamp	- •		
Auto Hold / Auto Stable / Real Time	• •		
Offset & Average Slope Display	• •		
Calibration Alarm (1 to 90 days)	• •		
Auto Shut-Off (1 to 30 mins.)	• •		
Electrode Status	• •		
Diagnostic Messages	• •		
PC Communication*1	-		
Printer Communication*2	-	•	
Meter Inputs	BNC, phono		
Display	Custom LCD with backlight		
Housing	IP67, shock & scratch resistant, non-slip		
Power Requirement	2 × AA batteries		
Battery Life	> 500 hours		
Dimensions	160 (L) × 80 (W) × 40.60 (H) mm		
Weight	Approx. 260 g (with batteries) / 216 g (without batteries)		

Specifications

- * Meter kit with NIST pH buffers is available upon request. Add 'N' suffix to order code.
- *1 Via PC (USB) cable PN 3200779639
- *2 Via Printer (RS232) cable PN 3200779638

HORIBA Advanced Techno

31, Miyanonishi-cho, Kisshoin Minami-ku, Kyoto 601-8306, Japan http://www.horiba-adt.jp

For any questions regarding this product, please contact your local agency, or inquire from the following website. http://global.horiba.com/contact_e/index.htm

