



POCKIT Central Automated PCR System

*Fast In-House
Molecular Testing*




Sample In – Result out!

POCKIT Central In-house PCR

Improving Veterinary Care

Specifications

Dimensions:	31cm x 48cm x 40cm
Net Weight:	21kg
Test Principle:	iiPCR
Storage Condition:	10-40°C 55 ± 25% RH
Operating Temperature:	15-35°C
Display:	Touchscreen GUI
Interface / Connectivity:	USB Port x3
Number of Samples:	1-8 reactions per run, per instrument
Sample Input Volume:	200µl per reaction
Sample Type:	Various
Sample Identification:	Bar code reader supplied (2D)
Run Time:	85 minutes
Results Storage:	10.000 onboard
QC:	+ve controls available for all tests
Sterilisation:	Integrated, automatic UV cycle
Power Supply:	100 120/200 240 V AC, 50/60 Hz, 2A

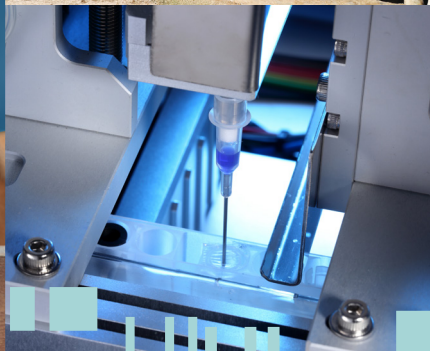

The **POCKIT Central** is a compact, fully integrated bench-top iiPCR instrument with a small footprint - suitable for fast in-house testing on any species or environment, improving diagnosis times and case management.

Set up and run your analysis in minutes.

POCKIT Central accommodates up to 8 samples of any combination, using a fully automatic walk-away system with proven iiPCR technology.

Run single tests or focussed panels including viruses, bacteria, parasites and protozoa.

The in-built touch screen and 2D scanner enable user-friendly input and operation, with direct display of test results and storage of 10,000 records.

All sample handling is fully automated, combining nucleic acid extraction, amplification and detection.

No additional equipment or specialist staff required!


Following each test run, the **POCKIT Central** initiates a self-cleaning cycle, using ultraviolet light to ensure it is ready to process the next batch of samples in quick succession.

The **POCKIT Central** in-house PCR analyser...


Sample in - Result Out

...within 85 minutes!

POCKIT Central

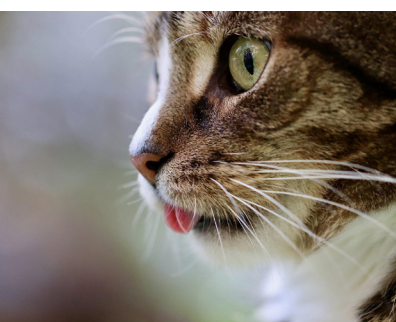

■ Preloaded ready-to-use cartridges

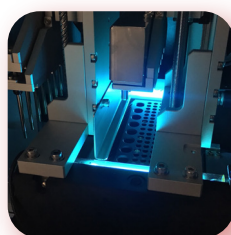

■ Built-in UV lamp for decontamination



■ Barcode scanning available



■ Result printing available

