



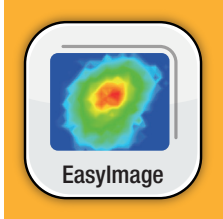

## EasyImage™: Makes Raman Imaging Easy

Simplicity taken  
to the next level




**EasyImage** is HORIBA's new imaging tool for Raman microscopy. Allowing quick and easy acquisition of meaningful chemical images of your samples without needing a high level of Raman expertise.

**EasyImage App** simplifies user's tasks. Meaning fewer user interventions and more time for productive work. You can now benefit from all of HORIBA's Raman know-how for your measurements!


**EasyImage App** offers an optimized analytical workflow for anyone, even busy users!

### Optimized Analytical Workflow


User friendly interface using **intelligent acquisition** and **auto-enhancing of Raman images!**

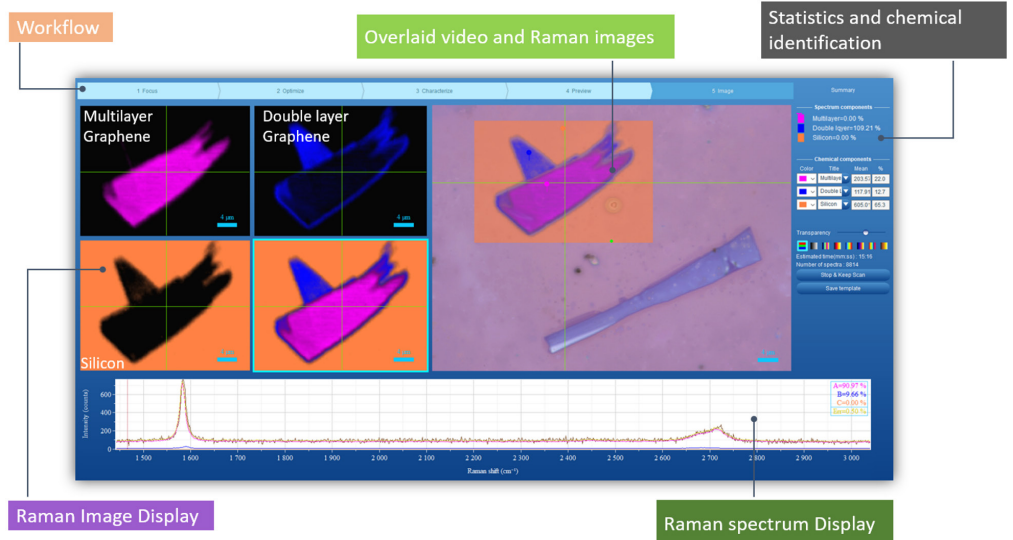


## Easy Image Workflow:

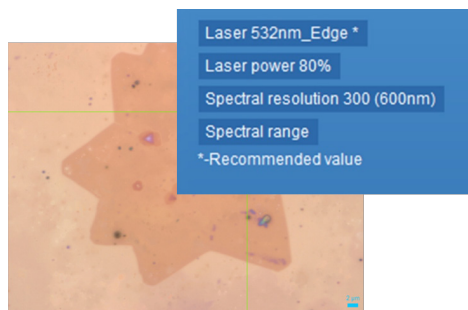


Automatic acquisition parameter selection (spectral range, laser intensity, contrast...)

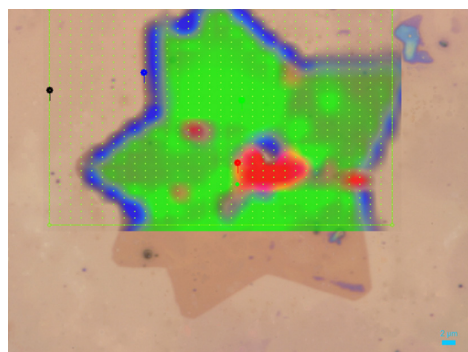

Mapping area selection

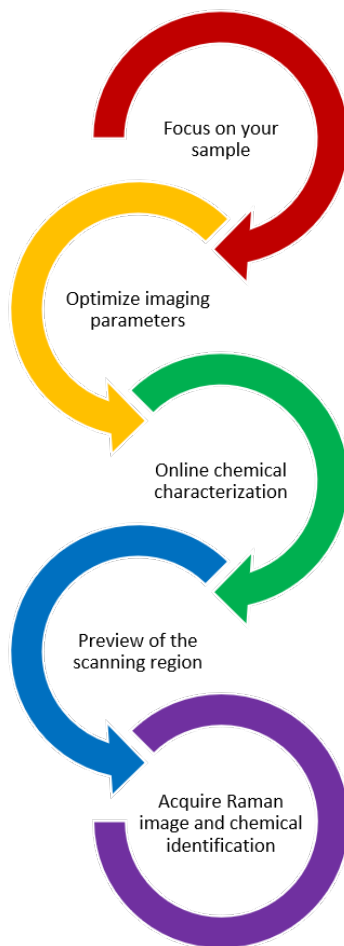

Targeting and focusing on X,Y

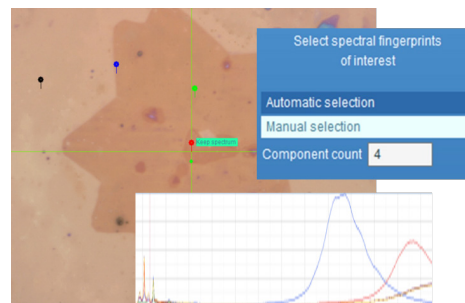

Chose the fingerprint!



Get your image!

## Easy Image main features:

- Automated objective detection for **LabRAM Soleil** only
- Intelligent spectral range selection
- Automated and manual spectral survey directly on the video image
- On-the-fly customized spatial resolution
- Automated background removal/using FLAT correction
- Automated spectral identification/with integrated KnowItAll® database
- Automated loading generation for chemometrics by *MVAPlus*
- Customizable templates and complete analysis report

**HORIBA**  
Scientific

[info.sci@horiba.com](mailto:info.sci@horiba.com)

**USA:** +1 732 494 8660  
**UK:** +44 (0)1604 542 500  
**China:** +86 (0)21 6289 6060  
**Taiwan:** +886 3 5600606

**France:** +33 (0)1 69 74 72 00  
**Italy:** +39 06 51 59 22 1  
**India:** +91 (80) 4127 3637  
**Brazil:** +55 (0)11 2923 5400

[www.horiba.com/scientific](http://www.horiba.com/scientific)

**Germany:** +49 (0) 6251 8475 0  
**Japan:** +81(75)313-8121  
**Singapore:** +65 (6) 745-8300  
**Other:** +33 (0)1 69 74 72 00

**20** YEARS  
1919-2019  
**HORIBA** JOBIN YVON