

LY-9610 Imaging Unit Accessory


See your particles!


Observe your particles in real time as they circulate in the Partica LA-960V2 wet system. Better understand your dispersions and extract image based shape and size information to complement your laser diffraction results. The small integrated unit fits neatly inside the Partica LA-960V2 without increasing the instrument footprint.

Example applications include:

- · Identify bubbles
- · Verify existence of oversize particles
- · Quantify particle shape

The separate software offers easy operation and full data export capabilities.

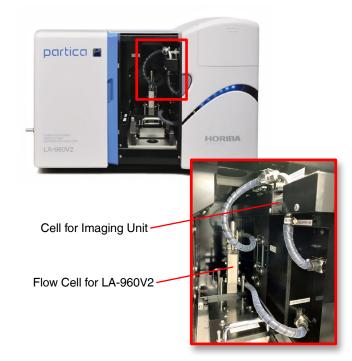


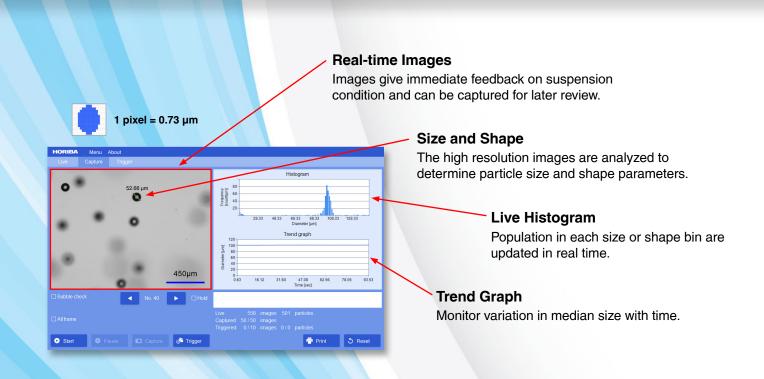

Features and Benefits

- Identify trace amounts of large particle contaminants or aggregates.
- Bubble detection.
- Improve your methods by developing an intuitive understanding of the effect of stirring, sonication, surfactant addition, etc.
- Detect large particles automatically.
- Shape distribution data (aspect ratio, circularity, length, width)*.
- Detailed analysis of captured images.
- Monitor trends in size and shape with trend graph and histogram.
- * Additional software required

Imaging Unit LY-9610

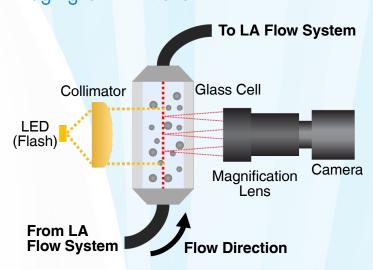

Specifications

Measureable Range for Image Analysis	9 - 1000 μm
Pixel Size	0.73 μm
Frame Rate	≥ 4 fps



https://bcove.video/2Bhuquw



The HORIBA Group adopts IMS (Integrated Management System) which integrates Quality Management System ISO9001, Environmental Management System ISO14001, and Occupational Health and Safety Management System OHSAS18001.

We have now integrated Business Continuity Management System ISO22301 in order to provide our products and services in a stable manner, even in emergencies



Please read the operation manual before using this product to assure safe and proper handling of the product.

- The specifications, appearance or other aspects of products in this catalog are subject to change without notice.
- Please contact us with enquiries concerning further details on the products in this catalog.
 The color of the actual products may differ from the color pictured in this catalog due to printing limitations.
 It is strictly forbidden to copy the content of this catalog in part or in full.
- The screen displays shown on products in this catalog have been inserted into the photographs through compositing.
 All brand names, product names, and service names in this catalog are trademarks or registered trademarks of their respective companies.

http://www.horiba.com

HORIBA, Ltd.	Japan
Head Office 2 Miyanohigashi-cho, Kisshoin, Minami-ku, Kyoto, 601-8510, Japan Phone: 81 (75) 313-8121 Fax: 81 (75) 321-5725	
HORIBA (China) Trading Co., Ltd.	China
Unit D, 1F, Building A, Synnex International Park, 1068 Tianshan Road, 200335, Shanghai, China Phone: 86 (21) 6289-6060 Fax: 86 (21) 6289-5553 Beijing Branch 12F, Metropolis Tower, No.2, Haidian Dong 3 Street, Be 100080, China Phone: 86 (10) 8567-9966 Fax: 86 (10) 8567-9066 Guangzhou Branch	
Room 1611 / 1612, Goldlion Digital Network Center, 138 Tiyu Road East, Guangzhou, 510620, China Phone: 86 (20) 3878-1883 Fax: 86 (20) 3878-1810	
HORIBA KOREA Ltd.	Korea
25, 94-Gif, Iljik-Ro, Manan-Gu, Anyang-Si, Gyeonggi-Do 13901, Korea Phone: 82 (31) 296-7911 Fax: 82 (31) 296-7913),

Singapore HORIBA Instruments (Singapore) Pte Ltd. 3 Changi Business Park Vista #01-01, Akzonobel House,

Singapore 486051 Phone: 65 (6) 745-8300 Fax: 65 (6) 745-8155 HORIBA Taiwan, Inc. Taiwan 8F.-8, No.38, Taiyuan St. Zhubei City, Hsinchu County 30265,

Phone: 886 (3) 560-0606 Fax: 886 (3) 560-0550

HORIBA India Private Limited

246, Okhla Industrial Estate, Phase 3 New Delhi-110020, India Phone: 91 (11) 4646-5000 Fax: 91 (11) 4646-5020 **Bangalore Office**

No.55, 12th Main, Behind BDA Complex, 6th sector, HSR Layout, Bangalore South, Bangalore-560102, India Phone: 91 (80) 4127-3637

HORIBA (Thailand) Limited Thailand

393, 395, 397, 399, 401, 403 Latya Road, Somdetchaopraya, Klongsan, Bangkok 10600, Thailand Phone: 66 (0) 2-861-5995 ext.123 Fax: 66 (0) 2-861-5200

East Office 850 / 7 Soi Lat Krabang 30 / 5, Lat Krabang Road, Lat Krabang, Bangkok 10520, Thailand

Phone: 66 (0) 2-734-4434 Fax: 66 (0) 2-734-4438

PT HORIBA Indonesia Indonesia

Jl. Jalur Sutera Blok 20A, No.16-17, Kel. Kunciran, Kec. Pinang Tangerang-15144, Indonesia Phone: 62 (21) 3044-8525 Fax: 62 (21) 3044-8521

HORIBA Vietnam Company Limited

Lot 3 and 4, 16 Floor, Detech Tower II, No.107 Nguyen Phong Sac Street, Dich Vong Hau Ward, Cau Giav District, Hanoi, Vietnam Phone: 84 (24) 3795-8552 Fax: 84 (24) 3795-8553

HORIBA Instruments Incorporated

9755 Research Drive, Irvine, CA 92618, U.S.A. Phone: 1 (949) 250-4811 Fax: 1 (949) 250-0924 HORIBA New Jersey Optical Spectroscopy Center 20 Knightsbridge Rd, Piscataway, NJ 08854, U.S.A. Phone: 1 (732) 494-8660 Fax: 1 (732) 549-5125

HORIBA Instruments Brasil, Ltda. Brazil

Rua Presbitero Plinio Alves de Souza, 645, Loteamento Multivias, Jardim Ermida II - Jundiai Sao Paulo - CEP 13.212-181 Brazil Phone: 55 (11) 2923-5400 Fax: 55 (11) 2923-5490

HORIBA FRANCE SAS France

16-18, rue du Canal, 91165, Longjumeau Cedex, France Phone: 33 (1) 69-74-72-00 Fax: 33 (1) 69-09-07-21

HORIBA Jobin Yvon GmbH Germany Neuhofstrasse 9, D_64625, Bensheim Phone: 49 (0) 62-51-84-750 Fax: 49 (0) 62-51-84-7520

HORIBA ITALIA Sri Italy

Via Luca Gaurico 209-00143, ROMA

Phone: 39 (6) 51-59-22-1 Fax: 39 (6) 51-96-43-34 UK **HORIBA UK Limited**

Kyoto Close Moulton Park Northampton NN3 6FL UK Phone: 44 (0) 1604 542500 Fax: 44 (0) 1604 542699

BJG052020

USA

