

SynapsePlus FIVS

High Speed Scientific CCD Camera

ELEMENTAL ANALYSIS
FLUORESCENCE
GRATINGS & OEM SPECTROMETERS
OPTICAL COMPONENTS
FORENSICS
PARTICLE CHARACTERIZATION
RAMAN
SPECTROSCOPIC ELLIPSOMETRY
SPR IMAGING

Front illuminated UV enhanced sensor,
-80°C (-95°C) 2048 x 512 pixels

Features and Benefits



- Deep thermoelectric cooling
- Fast acquisition rate of more than 2200 spectra/s
- Ideal for low light level detection without etaloning
- Excellent linearity
- e2V Scientific Grade 1 CCD
- Lifetime vacuum warranty
- USB 2.0 interface: 100% guaranteed data integrity
- HORIBA SynerJY or LabSpec acquisition and analysis software
- LabVIEW VI's and SDK available

Primary Applications

Primarily chosen for broad spectrum analysis such as photoluminescence, it is also well suited for studying fine spectral features on a broad spectral background.

- Absorption
- Raman
- Transmission
- Photoluminescence
- Reflectance
- Fast kinetics


QE Curve, SynapsePlus FIVS CCD

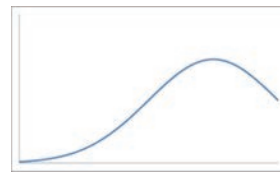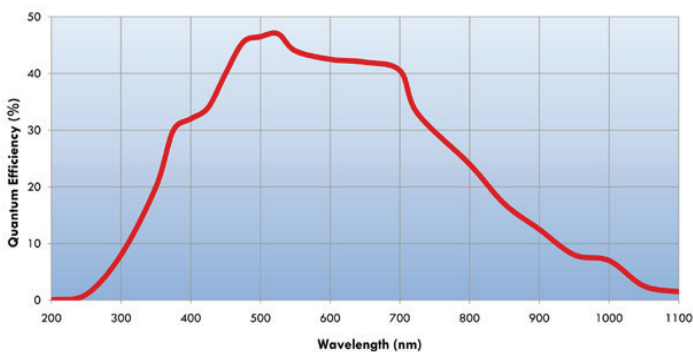

Specifications

Sensor	Front-Illuminated Visible (FIVS), Scientific Grade 1		
Active pixels	2048 x 512		
Pixel size	13.5 μm x 13.5 μm		
Image area	27.6 mm x 6.9 mm, 100% fill factor		
Pixel well depth	250,000 e ⁻		
Register well depth	1,000,000 e ⁻		
Non-linearity (Measured at all speeds per camera)	<0.40%		
Readout Noise (e⁻): Typ. (Max)	50 kHz: 2.5 (4)	1 MHz: 7 (12)	3 MHz: 11 (15)
Dark current at -80°C (e⁻/pixel/s)	<0.001		
Maximum spectra per second	Full vertical bin: 181 ROI mode 20 rows: 1795 ROI mode 8 rows: 2247		
Software-adjustable gain (e⁻/count)	Selectable from 1.0 to 13.5		
Digitization	16-bit ADC		
Vertical shift rates (μs)	9.3, 18.4, 36.5 software-selectable		
Cooling at +20°C	Air-cooled: -75°C (guaranteed), -80°C (typical) Liquid-cooling recirculator: -90°C (guaranteed), -95°C (typical)		
Power requirements AC-DC power supply (provided)	AC input 90-264 VAC, 47-63 Hz DC output, +9 V, 6.44 A maximum		
Physical dimensions	Length: 181 mm (7.14")	Width: 113 mm (4.45")	Height: 113 mm (4.45")
Physical weight	2.10 kg (4.64 lb)		

Mechanical Dimensions

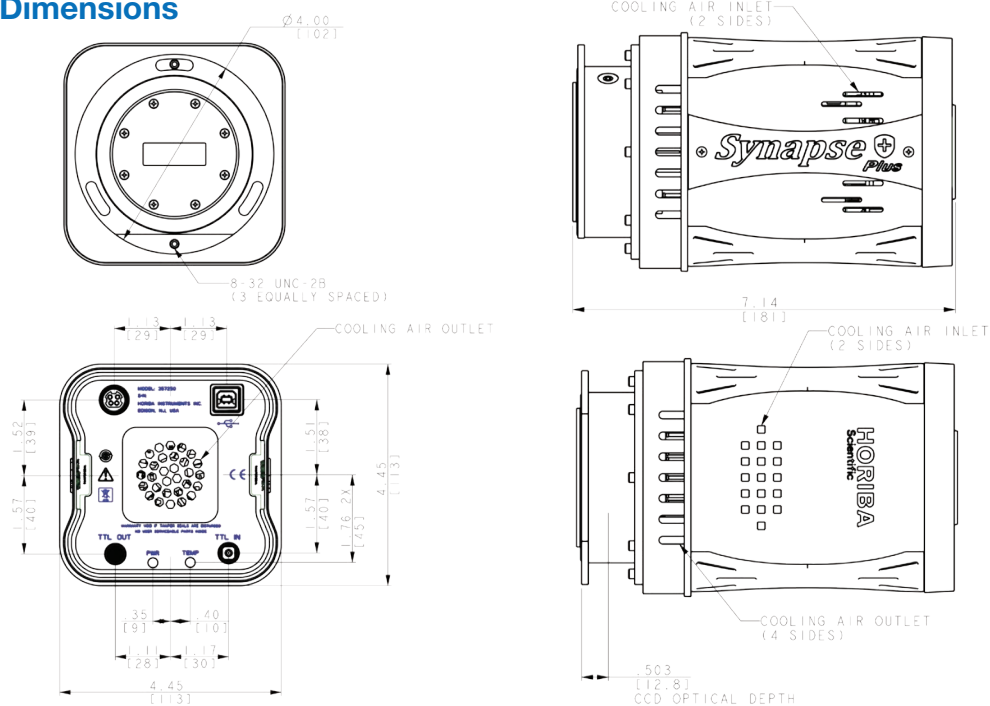

info.sci@horiba.com www.horiba.com/osd

USA: +1 732 494 8660
UK: +44 (0)20 8204 8142
China: +86 (0)21 6289 6060

France: +33 (0)1 69 74 72 00
Italy: +39 2 5760 3050
Brazil: +55 (0)11 2923 5400

Germany: +49 (0)6251 8475 0
Japan: +81 (0)3 6206 4721
Other: +1 732 494 8660

HORIBA
Scientific